

How often can I use my Philips Sonicare toothbrush per charge?

Do you want to know the operating time of your fully charged Philips Sonicare toothbrush? Find the operating time for your Philips Sonicare toothbrush.

Operating time for each toothbrush

Below operating times are based on a full charge and regular use (twice a day). It can take up to 24 hours before your Philips Sonicare toothbrush is fully charged. To avoid an empty battery, you can always leave your toothbrush on the charger between brushings.

Operating time from a full to empty

- DiamondClean Smart; up to 2 weeks
- DiamondClean; up to 3 weeks
- Flexcare Platinum Connected; up to 2 weeks
- Flexcare and Flexcare+; up to 3 weeks
- ProtectiveClean; up to 2 weeks
- 3 Series; up to 3 weeks
- 2 Series; up to 2 weeks
- Healthywhite+; up to 3 weeks
- Essence+; up to 10 days
- Essence; up to 2 weeks

Can I use my Philips Sonicare charger with other toothbrushes?

Philips Sonicare chargers are designed to fit the matching toothbrush handle, however, some models are interchangeable. Find out below if your charger is compatible with other Philips Sonicare toothbrushes.

Interchangeable Philips Sonicare chargers

Below you will find a list of Philips Sonicare toothbrushes that are compatible with other chargers:

- DailyClean toothbrush
- Elite toothbrush
- EasyClean toothbrush
- Essence+ toothbrush
- ExpertClean toothbrush
- ExpertResults toothbrush
- Flexcare Platinum Connected toothbrush
- FlexCare Platinum toothbrush
- FlexCare+ toothbrush
- FlexCare toothbrush
- HealthyWhite toothbrush
- HealthyWhite+ toothbrush
- ProtectiveClean toothbrush
- 2 Series Plaque Removal toothbrush
- 3 Series Gum Health toothbrush
- Sonicare for Kids
- Philips Sonicare AirFloss

The DiamondClean and DiamondClean Smart chargers

The DiamondClean and DiamondClean Smart chargers are not interchangeable with each other or with other Philips Sonicare toothbrush handles. Even if the charging light on the handle blinks, it will not charge properly. We advise to always use the charger that comes with your Philips Sonicare DiamondClean and DiamondClean Smart toothbrush.

Chargers can look similar. To verify you are using the correct charger, see the table below. Each handle is only compatible with the charger base below the handle code ID.

Handle Name	DiamondClean	DiamondClean Smart	DiamondClean 9000
Code on bottom of handle	HX93, HX934, HX935, HX936, HX937	HX992, HX993	HX991
Code on bottom of charger base	HX9100	Type: CBA1001, CBA2001, CBA3001	Type: CBA5001, CBA4001

Which devices are compatible with my Philips Sonicare App?

Your Philips Sonicare App supports every smartphone that runs on Android 6.0 and higher and every iPhone that runs on iOS 11.0 and higher. You can download the app by clicking here on [Google Play Store](#) and or the [iOS App Store](#).

What are the differences in modes of my Philips Sonicare Toothbrush?

Some Philips Sonicare Toothbrush models come with different brushing modes. You can find here an overview of what these modes mean. Please note that these modes may vary per toothbrush model.

For specific details regarding your particular toothbrush model refer to the user manual.

Brushing modes

Clean Mode

This mode has a higher frequency of vibrations in comparison to a standard mode for better cleaning. Brushing time is 2 minutes.

Deep Clean mode

You can choose this mode for extra cleaning. It is more powerful than any other mode. Brushing time is 3 minutes.

White mode

In this mode the total brushing time is two and a half minutes. The first 2 minutes are clean mode and the last half minute is polishing mode.

Polish mode

This mode lasts one minute. The toothbrush runs at a lower frequency of vibrations than in the standard cleaning mode.

Gum Care mode

This mode has a 1-minute massage cycle after the standard 2-minute clean mode to clean your

gums. Brushing time is 3 minutes.

Gum Health mode

This mode has an extended massage cycle after the standard 2-minute clean mode to clean your gums. Brushing time is 3:20 minutes

Sensitive mode

This mode has a lower frequency of vibrations in comparison to the standard mode. Brushing time is 2 minutes.

Tongue Care mode


This is a 20-second mode designed to clean your tongue, together with the Tongue Care brush head. Brushing time is 20 seconds

Changing the intensity settings

Your Sonicare Toothbrush gives you the option to choose between 3 different intensity levels: low, medium and high. Based on the brush head you attach, the intensity will be selected automatically. If you wish to change to a different intensity, you can do so by pressing the mode/intensity button while brushing.

Note: This feature works for Philips Sonicare DiamondClean Smart, ExpertClean, Expert Results, Protective Clean 5100 and 6100 Series.

For toothbrushes with an Intensity Level Panel on the handle, press the intensity panel up or down to adjust the intensity level. Check your user manual for more details.

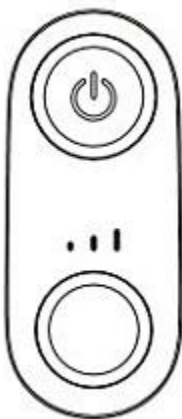

Which accessories are available for my Philips Sonicare Toothbrush?

If you are looking for new parts and accessories for your Philips Sonicare Toothbrush, please follow the steps below to find out which accessories are compatible with your Philips Sonicare Toothbrush.

Accessories for Philips Sonicare Toothbrush

Follow these steps to see which accessories are available for your Philips Sonicare Toothbrush:

1. Go to [Philips' parts and accessories page](#)
2. Use the search bar to search by model number or navigate to your product
3. Once on the electric toothbrush overview page, you can use the filters on the left to specify your search

Please contact us if you need further assistance.

Can I use my Philips Sonicare Toothbrush without the Sonicare app?

You can use your Philips Sonicare Toothbrush without the Sonicare App, but we advise you to use it as much as possible to maintain great oral health. When you connect your Philips Sonicare Toothbrush to the App, you will receive personalized guidance and tips to improve your brushing habits, leaving you with a healthy mouth and smile. Find out more about this in the following lines.

Real-time brushing guidance

The Philips Sonicare App has the guided brushing feature that will help you to improve and maintain your oral health by:

- * Tracking how often and how long you have brushed each time,
- * Highlighting missed spots on your teeth,
- * Notifying you to replace the brush head in time.

Note: The Philips Sonicare Flexcare Connected and DiamondClean Smart Toothbrush can store up to about 28 brushing sessions of data. To avoid losing your brushing data, regularly connect your toothbrush with the Sonicare App.

Can I exchange points in the Philips Sonicare App?

No, you cannot exchange points in the Philips Sonicare App at this time.