



The SafeWander Bed Exit Alarm Kit package contains:

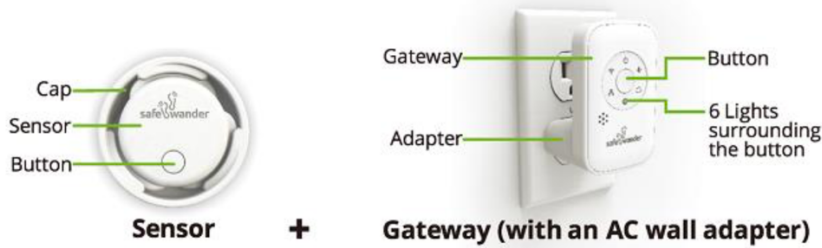

Suggested sensor locations

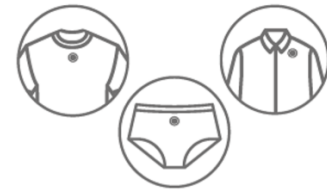

A QUICK START

Follow the steps below to create your user account and register the sensor and gateway.

Step I: Download the “SafeWander” app at the App Store or GooglePlay. SafeWander is compatible to Apple or Android mobile devices including smartphones or tablets.

Step II: One-Time Registration: Open the App, follow its instructions to

- 1) Create a user account using your email address and a password;
- 2) Add your loved one’s name (or nickname);
- 3) Register your Sensor;
- 4) Plug in the Gateway to an outlet ideally less than 10 ft away from your loved one’s bed. Link the Gateway to your WiFi (you need the password of your WiFi router.)

Caution: Choose 2.4 GHz WiFi. SafeWander is NOT compatible with 5.0 GHz WiFi.

Step III: Attach the Sensor to Clothing

Watch a short video on <http://www.safewander.com/support> to learn the easy and correct way to attach the Sensor to clothing.

Caution: SafeWander triggers an alert on your phone when the Sensor detects an angle change from a lying down to a sitting up position in the flexion direction. Turning left or right on the bed will not trigger an alert. So it is critical to make the Sensor move with the person’s body by *using a tight-fitting shirt or underwear* and to position the Sensor in such an orientation that *the SafeWander logo on the Sensor is approximately parallel with the shoulder line*.



Step IV: Test SafeWander® Alert:

- 1) On My Patients dashboard of the app, turn on the Alert Switch next to your loved one’s name.
- 2) When the Sensor is lying flat with the SafeWander logo facing up, the status should show “Sleeping” on the app.
- 3) Manually lift the Sensor from the lying to the sitting up position in the flexion direction (as if the SafeWander logo is the shoulder line), and you phone should beep.
- 4) Touch the screen to stop the beeping alert;
- 5) Tap “Confirm” (if it is a real event) or “Ignore” (if it is false)



Step V: Test SafeWander® in Two Different Wireless Modes:

SafeWander® is designed with dual wireless modes, and you may want to test both of them.

The Bluetooth Mode: When your phone (with its Bluetooth turned on) is in the same room as the Sensor, they will automatically link without routing via the Gateway. In this mode, the sensor light (at the 1-o'clock position) on the Gateway is not lit. When you travel with your loved one and stay in the same hotel room, you only need to take the Sensor, not the Gateway.

The WiFi Mode: If the Bluetooth of your phone is turned off or if your phone is moved away from the Sensor, the Sensor will automatically link to the nearby Gateway. In this mode, the sensor light on the Gateway is lit in green. These two wireless modes are automatically switched, but it may take a few seconds to switch. If you want to test the WiFi mode when you phone is close to the Sensor, you may turn off the Bluetooth on your phone.

USE OF SAFEWANDER®

Use as a Bed Alarm: After your loved one goes to bed, turn on the Alert Switch on your app to start monitoring. Make sure that the status on your app shows “Sleeping”

When you do not want to receive the beeping alert, simply turn off the Alert Button on the app. You can still monitor the status of the patient: Sleeping or Awake.

Use as a Chair Alarm: You can attach the Sensor to pants (on top of a thigh) to monitor chair exit.

Turn on/Turn off the Sensor: The power of the Sensor is turned on or off by pressing and holding the button on the Sensor for a few seconds until the light on the Sensor flashes.



Check if the Sensor is On or Off: Press the button on the Sensor once. If the light flashes, the Sensor is on. Otherwise, the Sensor is off.

It is recommended to keep the Sensor on even when you are not using it, unless you will not use it for 2 weeks or longer.

History Dashboard: Each bed-exit event is recorded. View the log in the History dashboard.

SIMULTANEOUS MONITORING BY MULTIPLE CAREGIVERS

Multiple caregivers can simultaneously monitor one sensor. *They must sign in to the same user account using the same password.* A sensor can only be registered to one user account.

SIMULTANEOUS MONITORING OF MULTIPLE SENSORS

A caregiver can simultaneously monitor multiple sensors, i.e., multiple people. Tap My Patients dashboard on the app, then tap “+” to add an additional person, and register his/her Sensor.

For any further questions, visit FAQ at <http://www.safewander.com/support> or contact us via email at support@safewander.com and phone at (800) 844-4248.