Flash*&GO*

Instructions For Use

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Silk'n Flash&Go™ is a powerful device. As such, it should be used with special attention to safety


!Keep away from your eyes

Applying a Silk'n Flash&Go™ pulse to the eyes can cause eye injury

Do not use Silk'n Flash&Go™ on tanned skin or after recent sun exposure

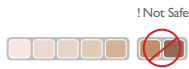Do not use Silk'n Flash&Go™ on naturally dark skin

Please read all warnings and safety precautions before use, and strictly follow them when using !™ Silk'n Flash&Go

Note! Silk'n Flash&Go™ may not be effective on naturally white and light hair

Quick Start

Plan your Silk'n Flash&Go hair removal for best results!

Treatments 1-4 – plan two weeks apart

Treatments 5-7 – plan four weeks apart

Treatments 8 + – treat as needed, until desired results are achieved.

Set up Silk'n Flash&Go!

1. Remove Silk'n Flash&Go™ HAIR REMOVAL UNIT and other components from box.
2. Verify that the DISPOSABLE LAMP CARTRIDGE is inserted properly at the tip of the device.
3. Plug the ADAPTOR cord into the Silk'n Flash&Go™ HAIR REMOVAL UNIT socket.
4. Plug the ADAPTOR into an electrical outlet. The STAND BY INDICATOR LIGHT will turn on. Your Silk'n Flash&Go™ is now ready to start.

Treat with Silk'n Flash&Go safely and properly!

5. Choose areas to treat that have not had recent exposure to the sun.
6. The skin of the treatment areas should be shaved, clean, dry and free of any powders, antiperspirants or deodorants. Do not wax, pluck or tweeze the hair out. Shaving is an important step to get the results you want!
7. Press the ON/STBY SWITCH. The POWER ON INDICATOR LIGHT (orange) will turn on and a fan sound (similar to the sound of a hairdryer) will start.
8. Approximately 3.5 seconds after pressing the ON/STBY SWITCH, the READY INDICATOR LIGHT (green) will turn on. The device is then ready for you to trigger the first pulse at the lowest energy level.
9. If this is your first treatment - Choose the energy level. Flash&Go offers 5 levels of energy from 1 (the lowest) to 5 (the highest). To choose the right energy level for your first treatment, test each treatment area separately:

- Place the **TREATMENT SURFACE** on the skin and press the trigger to deliver one pulse.
- If you did not experience abnormal discomfort press the + button once and deliver one pulse on a different spot at energy level 2.
- Continue testing at increased energy levels until you have reached the highest level that you feel comfortable with. If within one hour you experience no adverse effects, you can start your full treatment at that energy level.
- Repeat this test for each body part you intend to treat.

10. Set the energy level. If the energy level should be higher than the lowest, press the **ENERGY LEVEL SETTING BUTTONS** using the “-” or “+” to respectively decrease or increase the energy level, until the desired energy level is set and represented by 3 energy level indicator lights (For details see “Energy Level” chapter inside this user manual).

11. Apply the TREATMENT SURFACE to the skin, making sure the skin is spread evenly and smoothly. As soon as the TREATMENT SURFACE is in full contact with the skin, the READY INDICATOR LIGHT will start blinking.
12. Press the PULSE BUTTON. The device will first determine the color of your skin. If the color of the skin is light enough for safe application, the device will flash a pulse of light onto your skin, and the READY INDICATOR LIGHT will turn off. You will see a bright flash of light and simultaneously hear a subtle pop sound, which is a normal noise for the device. You will feel a mild sensation of warmth and tingling. Silk'n Flash&Go™ will immediately recharge for the next pulse. After 3.5 seconds the READY INDICATOR LIGHT will turn on again.
13. Remove the TREATMENT SURFACE from the treated area of skin.

14. Move the TREATMENT SURFACE to another area of skin. Use the pressure marks the device just made on your skin to guide you for proper positioning of the next pulse, avoiding both gaps and overlaps between pulses.

Try to avoid overlapping pulses!

Warning: Do not treat the same area of skin more than once per hair removal session!

Treating the same area of skin more than once per session increases the likelihood of adverse effects.

15. Repeat the process until the entire treatment area is covered.

**If your skin blisters or burns,
STOP USE IMMEDIATELY!**

Home Skinovations Ltd. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd. is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd.

Data is subject to change without notification.

Home Skinovations Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd.

Specifications are subject to change without notice.

Home Skinovations Ltd. Tavor Building, Shaar Yokneam POB 533, Yokneam 20692, ISRAEL www.silkn.com, info@silkn.com

Table of Contents

Silk'n Flash&Go Intended Use.....	p.7
Everything there is to know about Flash&Go.....	p.8
Silk'n Flash&Go device description.....	p.8
The disposable lamp cartridge.....	p.8
The skin color sensor.....	p.9
Energy levels.....	p.10
What to expect.....	p.11
Possible side effects.....	p.12
Plan your Flash&Go hair removal for best results!.....	p.15
Get started!.....	p.16
Treat with Silk'n Flash&Go safely and properly!.....	p.16
When Silk'n Flash&Go™ session has been completed.....	p.19
Skin care following hair removal session.....	p.20
To Replace a Disposable Lamp Cartridge.....	p.20
Safety with Silk'n Flash&Go.....	p.21
Troubleshooting.....	p.29
Customer Service.....	p.30
Specifications.....	p.31
Labeling.....	p.32
EC Declaration of Conformity.....	p.33

Before using Silk'n Flash&Go™ for the first time, please read this User Manual in its entirety.

Pay particular attention to sections on device use procedures, device operation, and after-use procedures.

We recommend you re-familiarize yourself with this User Manual before each use of Silk'n Flash&Go™

Silk'n Flash&Go™ is a powerful electrical device. As such, it should be used with special attention to safety.

Please read all warnings and safety precautions before use, and strictly follow them when using Silk'n Flash&Go™.

Silk'n Flash&Go Intended Use

- Silk'n Flash&Go is a light based device intended for the removal of unwanted hair.
- It was designed for safe 'at-home' use.
- It will effect long term, or permanent body hair reduction (permanent hair reduction is defined as a long-term stable reduction in the number of hairs re-growing after a treatment regime.)
- Ideal body areas for Silk'n Flash&Go use include the underarms, bikini line, arms, legs, face, back, shoulders and chest.
- The Silk'n Flash&Go is not suitable for naturally dark skin tones, tanned skin or after recent sun exposure.
- Silk'n Flash&Go is not as effective on blonde, white or light hair.

Everything there is to know about Flash&Go

Silk'n Flash&Go device description:

Your Silk'n Flash&Go consists of a **HAIR REMOVAL UNIT**(1) with **CONTROL PANEL** (2), **PULSE BUTTON** (3) and a **DISPOSABLE LAMP CARTRIDGE** (4).

Upon opening the Silk'n Flash&Go™ package, you will find the following parts:

- Silk'n Flash&Go™ Hair Removal Unit
- An Electrical Adaptor
- Extra Disposable Lamp Cartridge
- This User's Manual
- An Instructional DVD

Disposable Lamp Cartridge

The Silk'n Flash&Go™ DISPOSABLE LAMP CARTRIDGE can fire 1000 light pulses covering an area of 4 cm² with each pulse. Pulse intensity is determined only according to the energy level setting of the device. There is no decrease of energy during the usable lifetime of the LAMP CARTRIDGE.

When a **DISPOSABLE LAMP CARTRIDGE** has reached 90% of its lifetime the **REPLACE CARTRIDGE INDICATOR LIGHT** on the **CONTROL PANEL** will start blinking, indicating that the **DISPOSABLE LAMP CARTRIDGE** should soon be replaced.

As the **DISPOSABLE LAMP CARTRIDGE** usage is getting closer to 100%, the blink on-time of the **REPLACE CARTRIDGE INDICATOR LIGHT** will grow longer.

Once all 1000 pulses in a **DISPOSABLE LAMP CARTRIDGE** have been used the **REPLACE CARTRIDGE INDICATOR LIGHT** will be steadily on, indicating that further pulses cannot be emitted. In order to continue the hair removal session, it will be necessary to replace the **DISPOSABLE LAMP CARTRIDGE**.

For full instructions please see chapter: "To Replace a Disposable Lamp Cartridge".

Note: A **DISPOSABLE LAMP CARTRIDGE** should be replaced if large spots appear inside it or if its **TREATMENT SURFACE** is broken.

The Skin Color Sensor

- Light-based hair removal on darker skin tones can result in adverse effects such as burns, blisters, and skin color changes (hyper- or hypo-pigmentation).
- To prevent such misuse, a unique **SKIN COLOR SENSOR** in

Silk'n Flash&Go™ measures the treated skin complexion at the beginning of each session and occasionally during the session. If the SKIN COLOR SENSOR detects a skin tone that is too dark for Silk'n Flash&Go™ application, the device will automatically stop emitting pulses.

- If you see no light pulse and the SKIN TONE WARNING INDICATOR LIGHT blinks, this is an indication that your skin tone, measured by the SKIN COLOR SENSOR, is too dark for safe application. Try using the device on a different body part or contact the Silk'n Flash&Go™ support.

Energy level

The energy level determine the intensity of the Silk'n Flash&Go™ light pulse delivered to your skin, from the lowest level (-) to the highest level (+).

Flash&Go offers 5 levels of energy from 1 (the lowest) to 5 (the highest), represented by 3 energy level indicator lights:

Level 1 (lowest)	the left indicator light is solid green
Level 2	the left indicator light is solid green and the central indicator light is blinking green
Level 3	both left and central indicator lights are solid green

Level 4	the left and central indicator lights are solid green and the right indicator light is blinking green
Level 5 (highest)	all 3 indicator lights are solid green

When Silk'n Flash&Go™ is turned on its energy level will automatically be reset to the lowest energy level. Only one ENERGY LEVEL INDICATOR LIGHT will be on.

As energy level increases, so does hair removal results as well as the risk of possible side effects and complications. To achieve efficient hair removal from your first treatment, while avoiding side effects, you should test your skin tone before you start, and determine the optimal energy level for you. See instructions on how to test your skin tone on chapter: “Treat with Silk'n Flash&Go safely and properly!”.

What to expect when treating with the Silk'n Flash&Go

For many people, using Silk'n Flash&Go™ may be their first experience with a light-based device designed for home use. Silk'n Flash&Go™ is simple to use, and hair removal sessions are quick and easy. During a Silk'n Flash&Go™ session it is normal to experience and feel:

- **A Fan Noise** – The cooling fan in Silk'n Flash&Go™ makes noise similar to a hairdryer.
- **A Pop Sound with Each Pulse** – When a pulse of HPL™ light is activated, it is normal to hear a subtle pop sound simultaneously with the flash of light.
- **Moderate Pressure of the TREATMENT SURFACE** – This is necessary and helpful for placement of adjacent pulses of light, and is part of the unique safety feature of Silk'n Flash&Go™.
- **A Sensation of Warmth and Tingling** – During each pulse of light it is normal to feel a sensation of warmth and tingling from the light energy. You may feel the warmth for up to an hour after your Silk'n Flash&Go session.
- **Some Mild Red or Pink Color** – During and just after your Silk'n Flash&Go™ session it is not uncommon to see some very mild, pink-like color of the skin. This is usually most noticeable around the hairs themselves. However if you see full redness of the skin, blistering or burns stop use of Silk'n Flash&Go™ immediately.

Possible side effects using Silk'n Flash&Go

When used according to the instructions, side effects and complications associated with use of Silk'n Flash&Go™ are uncommon. However every cosmetic procedure, including those designed for home use, involves some degree of risk. Therefore it is important that you understand and accept the risks and complications that can occur with pulsed light hair removal systems designed for home use.

Adverse Reaction	Degree of Risk
Stacking or Overlapping of multiple pulses on the same skin spot	Minor
Minor Skin Discomfort – Although home pulsed light hair removal is generally very well-tolerated, most users do feel some mild discomfort during use, usually described as being a mild stinging sensation on the treated skin areas. The stinging sensation usually lasts during the time of the application itself or for a few minutes thereafter. Anything beyond this minor discomfort is abnormal and means that either you should not continue to use Silk'n Flash&Go™ because you are unable to tolerate the hair removal application, or that the energy level setting is too high.	Minor

Adverse Reaction	Degree of Risk
<p>Skin Redness – Your skin may become red right after using Silk'n Flash&Go™ or within 24 hours of using Silk'n Flash&Go™. Redness generally clears up within 24 hours. See your doctor if redness does not go away within 2 to 3 days.</p>	<p>Minor</p>
<p>Increased Sensitivity to Skin – The skin of the treated area is more sensitive so you may encounter dryness or flaking of the skin.</p>	<p>Minor</p>
<p>Pigment Changes – Silk'n Flash&Go™ targets the hair shaft, in particular the pigmented cells in the hair follicle and the hair follicle itself. Nevertheless there is risk of temporary hyperpigmentation (increased pigment or brown discoloration) or hypopigmentation (whitening) to the surrounding skin. This risk of changes in skin pigmentation is higher for people with darker skin tones. Usually discoloration or changes to skin pigment are temporary and permanent hyperpigmentation or hypopigmentation rarely occur.</p>	<p>Rare</p>

Adverse Reaction	Degree of Risk
<p>Excessive Redness and Swelling – In rare cases treated skin may become very red and swollen. This is more common in sensitive areas of the body. The redness and swelling should subside within 2 to 7 days and should be treated with frequent applications of ice. Gentle cleansing is OK, but one should avoid exposure to sun.</p>	Rare

Plan your Flash&Go hair removal for best results!

Typical hair removal plan*:

- Hair removal sessions 1-4 – schedule treatments two weeks apart.
- Hair removal sessions 5-7 – schedule treatments four weeks apart.
- Hair removal sessions 8+ – treat as needed until long term results are achieved**.

* This is the recommended treatment schedule that has proven to produce the best results, but you may plan your personal treatment schedule differently and still get satisfactory results. Results vary from person to person therefore treatment schedules may also vary.

**A typical full hair growth cycle may take 18 months. During this time multiple Flash&Go sessions may be required in order to achieve permanent hair removal.

** Due to hormonal or other physiological changes dormant hair follicles may become active. Maintenance hair removal sessions with Silk'n Flash&Go™ may be required from time to time.

Get started!

1. Remove Silk'n Flash&Go™ HAIR REMOVAL UNIT and other components from box.
2. Verify that the **DISPOSABLE LAMP CARTRIDGE** is inserted properly at the tip of the device and that it is not damaged.
3. Plug the ADAPTOR cord into the Silk'n Flash&Go™ HAIR REMOVAL UNIT socket.
4. Plug the ADAPTOR into an electrical outlet. The **STAND BY INDICATOR LIGHT** will turn on.

Your Silk'n Flash&Go™ is now ready to start.

Treat with Silk'n Flash&Go safely and properly!

Choose areas to treat that have not had recent exposure to the sun.

5. The skin of the treatment areas should be shaved, clean, dry and free of any powders, antiperspirants or deodorants. Do not wax, pluck or tweeze the hair out. Shaving is an important step to get the results you want!
6. Press the ON/STBY SWITCH. The POWER ON INDICATOR LIGHT

(orange) will turn on and a fan sound (similar to the sound of a hairdryer) will start.

7. Approximately 3.5 seconds after pressing the ON/ STBY SWITCH, the READY INDICATOR LIGHT (green) will turn on. The device is then ready for you to trigger the first pulse at the lowest energy level.
8. If this is your first treatment - Choose the energy level. Flash&Go offers 5 levels of energy from 1 (the lowest) to 5 (the highest). To choose the right energy level for your first treatment, test each treatment area separately:

- Place the TREATMENT SURFACE on the skin and press the trigger to deliver one pulse.
- If you did not experience abnormal discomfort press the + button once and deliver one pulse on a different spot at energy level 2.
- Continue testing at increased energy levels until you have reached the highest level that you feel comfortable with. If within one hour you experience no adverse effects, you can start your full treatment at that energy level.
- Repeat this test for each body part you intend to treat.

9. Set the energy level. If the energy level should be higher than the lowest, press the ENERGY LEVEL SETTING BUTTONS using the

“-“ or “+” to respectively decrease or increase the energy level, until the desired energy level is set and represented by 3 energy level indicator lights (For details see “Energy Level” chapter inside this user manual).

10. Apply the TREATMENT SURFACE to the skin, making sure the skin is spread evenly and smoothly.
As soon as the TREATMENT SURFACE is in full contact with the skin, the READY INDICATOR LIGHT will start blinking.
11. Press the PULSE BUTTON. The device will first determine the color of your skin.
If the color of the skin is light enough for safe application, the device will flash a pulse of light onto your skin, and the READY INDICATOR LIGHT will turn off.
You will see a bright flash of light and simultaneously hear a subtle pop sound, which is a normal noise for the device. You will feel a mild sensation of warmth and tingling.
Silk'n Flash&Go™ will immediately recharge for the next pulse. After 3.5 seconds the READY INDICATOR LIGHT will turn on again.
12. Remove the TREATMENT SURFACE from the treated area of skin.
13. Move the TREATMENT SURFACE to another area of skin. Use the pressure marks the device just made on your skin to guide you for proper positioning of the next pulse, avoiding both gaps and overlaps between pulses.

Silk'n Flash&Go™ pulses should be administered in rows, starting at one end of each row and progressing sequentially towards the other end. This technique allows better control of skin coverage, and helps you avoid treating the same area more than once or overlapping skin areas.

When applied to the skin, the Silk'n Flash&Go™ TREATMENT SURFACE is designed to create temporary pressure marks on the treated area. These visible marks can be used for proper positioning of the next pulse.

Try to avoid overlapping pulses!

Warning: Do not treat the same area of skin more than once per hair removal session!

Treating the same area of skin more than once per session increases the likelihood of adverse effects.

14. Repeat the process until the entire treatment area is covered.

If your skin blisters or burns, STOP USE IMMEDIATELY!

When Silk'n Flash&Go™ session has been completed

Turn Silk'n Flash&Go™ off by pressing the ON/STBY SWITCH. All the lit indicators will turn off and the STBY INDICATOR LIGHT will turn on (Be sure to remember the last energy level setting you used, as it will not be restored when turning Silk'n Flash&Go™ on again.) Unplug the ADAPTOR cord from the electrical outlet.

After each hair removal session it is recommended that you clean your Silk'n™ Flash&Go device, especially the TREATMENT SURFACE and store in its original box.

Use a dry, clean cloth and a specially formulated cleaner for electronic equipment to gently wipe Silk'n Flash&Go™ surface, especially the TREATMENT SURFACE. Never immerse Silk'n Flash&Go™ or any of its parts in water!

Skin care following hair removal session

Do not expose treated areas of skin to the sun. Be sure to carefully protect the treated skin with sunscreen, throughout the hair removal period and for at least 2 weeks following the last Silk'n Flash&Go™ session.

Note! After replacing the disposable lamp cartridge with a new one, always set the energy to be one level lower than your normal setting. Raise the energy level back, only if you do not experience abnormal discomfort during or after the treatment.

To Replace a Disposable Lamp Cartridge

- Press the ON/STBY SWITCH to turn Silk'n Flash&Go™ to STBY mode.
- Unplug the power adaptor from the electrical outlet.
- Grasp the used DISPOSABLE LAMP CARTRIDGE on both sides, pull it out of the socket and discard with normal trash.
- Unwrap a new DISPOSABLE LAMP CARTRIDGE.
- Push the new DISPOSABLE LAMP CARTRIDGE gently into place.

Safety with Silk'n Flash&Go

Silk'n Flash&Go is not designed for everyone. Certain conditions may prevent the use of Silk'n Flash&Go. Please read the following precautions and contraindications prior to use. If you are unsure the Silk'n Flash&Go is safe for you, please contact your local Silk'n domain or consult with your doctor or dermatologist.

Warning	Possible results	Safety Tip
<p>Do not use Silk'n Flash&Go on naturally dark skin complexion.</p> <p>Not Safe!</p>  <p>See color skin chart at the inner cover</p>	<p>Treating dark skin with Silk'n Flash&Go can result in adverse effects such as redness and skin discomfort</p>	<p>A unique SKIN COLOR SENSOR is embedded in Silk'n Flash&Go™ to measure the treated skin complexion at the beginning of each session and occasionally during the session. The SKIN COLOR SENSOR ensures that pulses will only be emitted on suitable skin tones.</p>

Warning	Possible results	Safety Tip
<p>Do not use Silk'n Flash&Go on tanned skin or after recent sun exposure (including artificial tanning machines).</p> <p>Do not expose treated areas to the sun.</p>	<p>Treating with Silk'n Flash&Go before or after any sun exposure can result in adverse effects such as redness and skin discomfort. Tanned skin particularly following sun exposure, contains large quantities of the pigment Melanin. This applies to all skin types and complexions, including those which don't seem to tan quickly. The presence of large quantities of Melanin exposes the skin to higher risk when using any method of light based hair removal.</p>	<p>Avoid exposure to the sun for 4 weeks before your Silk'n Flash&Go™ treatment!</p> <p>Sun exposure includes unprotected exposure to direct sunlight of over 15 minutes constantly, or unprotected exposure to diffused sunlight of over 1 hour constantly. To protect recently treated skin when exposed to sunlight, be sure to thoroughly apply sunscreen SPF 30 or higher, for 2 weeks after each hair removal session.</p>

Warning	Possible results	Safety Tip
<p>Never use Silk'n Flash&Go around or near the eyes.</p>	<p>Although, Silk'n Flash&Go™ can be used for the treatment of facial hair; extra caution should be used on the face to avoid the eyes.</p>	<p>Silk'n Flash&Go™ has a built-in safety feature for eye protection. It has been designed so that a light pulse can not be emitted when the TREATMENT SURFACE is facing open air. The safety switch enables treatment only when the TREATMENT SURFACE is in contact with the tissue.</p>
<p>Do not treat on tattoos or permanent makeup, dark brown or black spots, (such as large freckles, birth marks, moles or warts), nipples, genitals or lips.</p>	<p>Treating any area with the Silk'n Flash &Go that is dark in colour or has more pigment can result in adverse effects such as burns, blisters, and skin color changes (hyper- or hypo-pigmentation).</p>	<p>Cover area with material that will not absorb the light such as a white cloth or white medical tape.</p>

Warning	Possible results	Safety Tip
<p>Do not treat on active eczema, psoriasis, lesions, open wounds or infections (cold sores), abnormal skin conditions caused by systemic or metabolic diseases (diabetes for example)</p>	<p>Light based treatments can cause adverse reactions on already sensitive areas.</p>	<p>Wait for the effected area to heal before using Silk'n Flash&Go™.</p>
<p>A treatment area that has a history of herpes outbreaks.</p>	<p>Light based treatments can cause adverse reactions on already sensitive areas.</p>	<p>Consult your physician to received preventative treatment before using Silk'n Flash&Go™.</p>
<p>Do not treat with Silk'n Flash & go if you are pregnant or nursing.</p>	<p>Hormonal changes during pregnancy could increase skin sensitivity and the risk of adverse effects.</p>	
<p>Do not treat on areas where you may want hair to grow.</p>	<p>Results are permanent.</p>	

Warning	Possible results	Safety Tip
<p>Keep Silk'n Flash&Go away from water! Do not place or store Silk'n Flash&Go™ where it can fall, be pushed or placed into a tub, sink or any other vessel containing water.</p>	<p>This may cause severe electrocution.</p>	<p>Do not use Silk'n Flash&Go™: while bathing or if it becomes damp or wet. Unplug immediately if it has fallen into water.</p>
<p>Do not attempt to open or repair your Silk'n Flash&Go™ device.</p>	<p>Opening Silk'n Flash&Go™ may expose you to dangerous electrical components and to pulsed light energy, either of which may cause serious bodily damage and/or permanent eye injury. It may also damage your device and will void your warranty.</p>	<p>Please contact Silk'n™ Customer Service if you have a broken or damaged device in need of repair.</p>

Do not use Silk'n Flash&Go if:

If the DISPOSABLE LAMP CARTRIDGE is damaged, or if its TREATMENT SURFACE is cracked, chipped or missing.		Stop use immediately and contact Silk'n™ Customer Service.
the attachments or accessories are not recommended by Home Skinovations Ltd.		Stop use immediately and contact Silk'n™ Customer Service.
it has a damaged cord or plug.	Keep the power cord away from heated surfaces.	Stop use immediately and contact Silk'n™ Customer Service.
you see or smell smoke when it is in use.		Stop use immediately and contact Silk'n™ Customer Service.
it is not working properly or if it appears damaged.		Stop use immediately and contact Silk'n™ Customer Service.

Do not use Silk'n Flash&Go if:

its SKIN COLOR SENSOR is cracked, or broken.		Stop use immediately and contact Silk'n™ Customer Service.
the outer shell is cracked or is coming apart.		Stop use immediately and contact Silk'n™ Customer Service.

Avoid using Silk'n Flash&Go if any of the following applies:

History of keloidal scar formation.

If you suffer from epilepsy.

If you have a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.

If you have a history of skin cancer or areas of potential skin malignancies.

If you have received radiation therapy or chemotherapy treatments within the past 3 months.

If you have any other condition which in your physician's opinion would make it unsafe for you to be treated.

If your skin was treated with a physician-prescribed drug within the past 6 months, please consult with your physician.

Consult with your physician prior to using the Silk'n Flash&Go.

Troubleshooting

Do not attempt to open or repair your Silk'n Flash&Go™ device. Only authorized Silk'n™ repair centers are permitted to perform repairs.

<p>“My Silk'n Flash&Go™ does not start.”</p>	<ul style="list-style-type: none">• Make sure the POWER ADAPTOR is properly connected to the Flash&Go™ device.• Make sure the POWER ADAPTOR is plugged into an electrical outlet on the wall.• If the problem persists, contact your local Silk'n™ Customer Service Center.
<p>“A light pulse is not emitted when I press the PULSE BUTTON”</p>	<ul style="list-style-type: none">• Check the REPLACE CARTRIDGE INDICATOR LIGHT. If it is on, disconnect Silk'n Flash&Go™ and replace the DISPOSABLE LAMP CARTRIDGE.• Make sure that the READY INDICATOR is on.• If within 10 seconds the READY INDICATOR remains off turn Silk'n Flash&Go™ off and back on by pressing the ON/ STBY SWITCH twice.• Check the SYSTEM WARNING INDICATOR LIGHT.• If the light is on turn Silk'n Flash&Go™ off and back on by pressing the ON/STBY SWITCH twice.• If the problem persists, contact your local Silk'n™ Customer Service Center.

“A light pulse is not emitted when I press the PULSE BUTTON”

- Check the REPLACE CARTRIDGE INDICATOR LIGHT. If it is on, disconnect Silk'n Flash&Go™ and replace the DISPOSABLE LAMP CARTRIDGE.
- Make sure that the READY INDICATOR is on.
- If within 10 seconds the READY INDICATOR remains off turn Silk'n Flash&Go™ off and back on by pressing the ON/STBY SWITCH twice.
- Check the SYSTEM WARNING INDICATOR LIGHT.
- If the light is on turn Silk'n Flash&Go™ off and back on by pressing the ON/STBY SWITCH twice.
- If the problem persists, contact your local Silk'n™ Customer Service Center.

Customer Service

For more information about Silk'n™ products please enter your local Silk'n™ website, www.silkn.eu or www.silkn.com.

If your device is broken, damaged, in need of repair, or for any other user assistance, please contact your local Silk'n™ Customer Service:


For USA and Canada customers: 1-877-DO-SILKN / 1-877-367-4556,
contact@silkn.com

For customers in other countries: info@silkn.com

Specifications

Specs, Spot Size		1.3cm x 3cm [4cm ²]
Speed		1 pulse every 3.5 second: 1.7 cm ² /sec
Technology		Home Pulsed Light™
Max Energy Level		Max 5J/cm ²
Wavelength		475-1200nm
Charging Time/Power Source		Continuous operation
Electrical Requirements		100-240VAC, 2A
Time Needed to Treat Lower Legs		30 minutes
Operation and Safety		New Skin Color Sensor ensures use only on suitable skin complexions. Safety sensor tip enables maximum safety & control. Easy to use.
Package Size		22cm x 22cm x 13cm
Temperature	Operating	10 to 40°C
	Storage	-40 to +70°C
Relative Humidity	Operating	30 to 75%rH
	Storage	10 to 90%rH
Atmospheric pressure	Operating	700 to 1060hPa
	Storage	500 to 1060hPa

Labeling

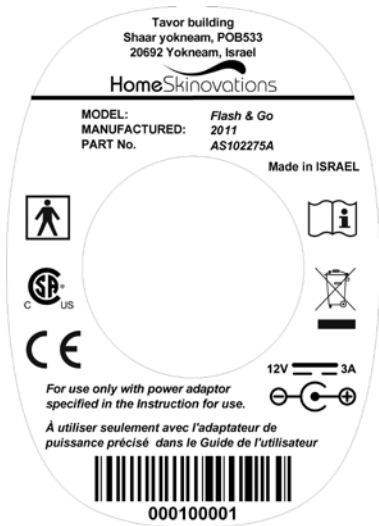Degree of protection against electric shock: type BF applied part

WEEE - Waste Electrical and Electronic Equipment

CE Mark

Follow operating instructions

CSA Mark for USA and Canada.

Degree of protection against ingress water: ordinary

This device is not suitable for use in the presence of flammable anesthetic mixture with air or with Oxygen or Nitrous Oxide.

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