

FULL BODY

LumaRx™

IPL SKIN BEAUTY SYSTEM



FULL-BODY HAIR REMOVAL SYSTEM

USE AND CARE GUIDE · IPL6800USA

REGISTER YOUR NEW PRODUCT TODAY!

By registering your new LumaRx™ product, you will receive these important benefits:

Proof of Ownership · Product Protection · Product Registration · Special Offers

Register online at: www.Luma-Rx.com

CONTENTS

PAGE

- 4 **Introducing the LumaRx™ Full-Body Hair Removal System**
- 5 **How does the LumaRx™ Full-Body Hair Removal System work?**
- 7 **Getting to know your LumaRx™ Full-Body Hair Removal System**
- 8 **IMPORTANT SAFETY INSTRUCTIONS**
- 9 **Servicing Of Double-Insulated Appliances**
- 10 **Facial cap**
- 11 **Changing treatment caps**
- 12 **Quick reference guide**
- 14 **Usage of the LumaRx™ Full-Body Hair Removal System (Indications)**
- 18 **Who shouldn't use the LumaRx™ Full-Body Hair Removal System (Contraindications)**
- 19 **WARNINGS**
- 20 **PRECAUTIONS**
- 23 **Risks**
- 25 **What to expect**

PAGE

- 27 **Operating Instructions**
 - Step 1: Prepare skin
 - Step 2: Turn on and unlock the LumaRx Full-Body Hair Removal System
 - Step 3: Select energy level
 - Step 4: Test a patch of skin
 - Step 5: Treat your skin
 - Other treatment tips
- 39 **User-Maintenance Instructions**
 - Cleaning the LumaRx Full-Body Hair Removal System
 - Replacing the Lamp Cartridge
 - Disposal of the LumaRx Full-Body Hair Removal System
- 44 **Troubleshooting**
- 51 **Contact information, specifications**
- 52 **Customer service, performance guarantee**
- 54 **Treatment cards**

INTRODUCING THE LUMARX™ FULL-BODY HAIR REMOVAL SYSTEM

Congratulations on the purchase of your new LumaRx Full-Body Hair Removal System! The LumaRx Full-Body Hair Removal System is a new breakthrough product that uses Intense Pulsed Light (“IPL”) to remove unwanted hair at-home. The LumaRx Full-Body Hair Removal System utilizes the same clinically proven technology that professional dermatologists and salons use, but optimized for personal treatment of smaller areas in the privacy and comfort of your own home.

Unlike most at-home hair removal methods, the LumaRx Full-Body Hair Removal System doesn’t cut or pull out hair; it actually uses light energy to work below the skin’s surface to disable active hairs at the root. You should gradually notice fewer hairs in the treatment area within a few weeks, with results generally improving over time.

The LumaRx Full-Body Hair Removal System is intended for use on:

- Female Facial Hair (cheeks, sideburns, chin, neck areas, and above the lips); with the Facial Cap.
- Male and Female Body Hair (underarms, bikini line, legs, arms, stomach and chest) with the Body Cap.

Note: The LumaRx Full-Body Hair Removal System isn’t for everyone. Please review the Contraindications, WARNINGS and PRECAUTIONS starting on page 19 to determine if the LumaRx Full-Body Hair Removal System is OK for you to use.

HOW DOES THE LUMARX™ FULL-BODY HAIR REMOVAL SYSTEM WORK?

Light-based hair removal is a well-established technology that has been in use worldwide for over 15 years. The LumaRx Full-Body Hair Removal System uses this technology to flash your skin with an intense pulse of light. The energy in the pulse works under the surface of your skin to target hairs at the root. Because the energy in the pulse of light targets dark coloring, it works best on darker hair. To avoid targeting (and possibly damaging) dark skin, you should use the LumaRx Full-Body Hair Removal System only on light-to-medium skin tones.

How Many Treatments are Needed?

Individual hairs go through a continuous growth cycle with phases of active growth and phases of rest. Like other laser or light-based treatments, the LumaRx Full-Body Hair Removal System affects only hairs that are in the active growth phase. Your initial regimen should be spread out over 1 month to ensure that you target the most number of active hairs (including those that recently cycled from the resting phase).

We recommend that you use the LumaRx Full-Body Hair Removal System once every other week over 1 month (for a total of 3 treatments) for each area.


If you have redness, blistering or other skin irritation in the treatment area, DO NOT perform any additional treatment. Wait until the skin returns to normal before using the LumaRx Full-Body Hair Removal System again. If the irritation does not go away within 3 days, contact Customer Service on page 51.

Use the **Treatment Record Cards** on pages 56-58 to help keep track of your treatments.

Results aren’t immediate, but you should have noticeably fewer hairs within 1 to 2 weeks after you start treatment. During this time, it may look like the hair is still growing, but it is actually being “pushed out.” You may or may not see these hairs fall out.

LUMARX™ FULL-BODY HAIR REMOVAL SYSTEM DIAGRAM REFERENCE

Figure 1

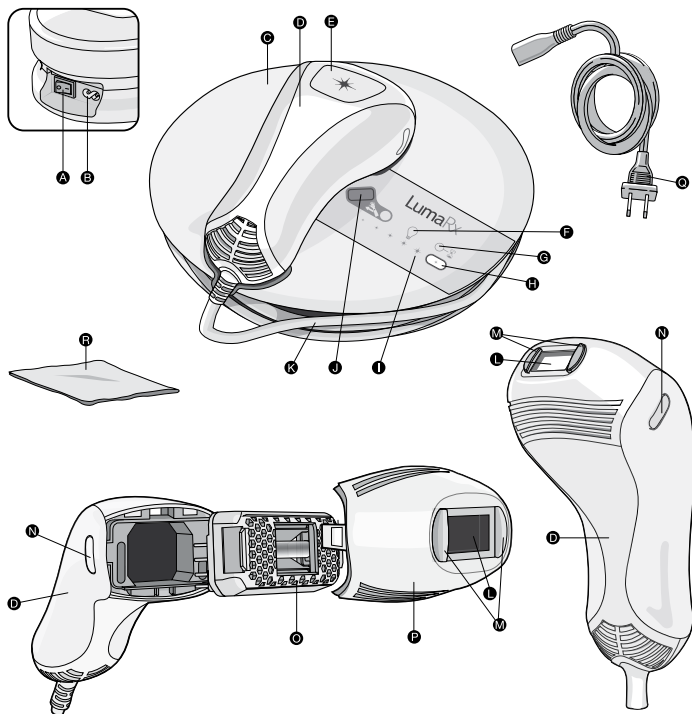

- Ⓐ Power Switch
- Ⓑ Power Socket
- Ⓒ Base Unit
- Ⓓ Hand Piece
- Ⓔ Flash Button
- Ⓕ Cartridge Status Indicator Display
- Ⓖ Cool Down Mode Indicator
- Ⓗ Energy Level Selection Button
- Ⓘ Energy Level Selection Display
- Ⓝ Skin Tone Tester
- Ⓚ Hand Piece Cord
- Ⓛ Flash Window
- Ⓜ Skin Contact Sensors
- Ⓝ Treatment Cap Release Buttons
- Ⓞ Replacement Cartridge
- Ⓟ Body Cap
- Ⓠ Power Cord
- Ⓡ Lint Free Cloth
- Ⓢ Facial Cap (not shown)

Getting to know your LumaRx™ Full-Body Hair Removal System

Inside your LumaRx Full-Body Hair Removal System box, you will find the BASE UNIT, POWER CORD, CLEANING CLOTH, and this USER'S MANUAL. The LumaRx Full-Body Hair Removal System comes with 2 TREATMENT CAPS (for Body and Facial areas) and a LAMP CARTRIDGE pre-assembled to the HAND PIECE. Each lamp cartridge provides 65,000 flashes before it needs to be replaced.

On the BASE UNIT, the control panel includes a skin tone sensor to test for appropriate skin tone before usage, and allows you to control the energy level that's effective and comfortable for you (see instructions for "Testing the LumaRx Full-Body Hair Removal System on a patch of skin" on pages 32-33). The Power Button is located next to the POWER SOCKET.

IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS BEFORE USING

KEEP AWAY FROM WATER

DANGER – To reduce the risk of electrocution:

1. Always unplug this appliance immediately after using.
2. Do not use while bathing.
3. Do not place or store appliance where it can fall or be pulled into a tub or sink.
4. Do not place in or drop into water or other liquid.
5. Do not reach for an appliance that has fallen into water. Unplug immediately.

WARNING – To reduce the risk of burns, electrocution, fire, or injury to persons:

1. An appliance should never be left unattended when plugged in.
2. Close supervision is necessary when this appliance is used by, on, or near

children or individuals with certain disabilities.

3. Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, or if it has been dropped, damaged, or dropped into water. Return the appliance to a service center for examination and repair.
5. Keep the cord away from heated surfaces. Do not wrap the cord around the appliance.
6. Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.
7. Never use while sleeping or drowsy.
8. Never drop or insert any object into any openings or hose.

9. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

10. Do not pull or carry by cord, use cord as a handle, close a door on cord, or pull cord around sharp edges or corners. Do not run Appliance over cord. Keep cord away from heated surfaces.

11. Do not unplug by pulling on cord. To unplug, grasp the plug, not the cord.

12. Do not handle plug or Appliance with wet hands.

13. Turn off all controls before unplugging.

SAVE THESE INSTRUCTIONS

SERVICING OF DOUBLE-INSULATED APPLIANCES

In a double-insulated appliance, two systems of insulation are provided instead of grounding. No grounding means is provided on a double-insulated appliance, nor should a means for grounding be added to the appliance. Servicing a double-insulated appliance requires extreme care and knowledge of the system, and should be done only by

qualified service personnel. Replacement parts for a double-insulated appliance must be identical to the parts they replace. A double-insulated appliance is marked with the words “DOUBLE INSULATION” or “DOUBLE INSULATED.” The symbol (square within a square) may also be marked on the appliance.


FACIAL CAP

Use to treat facial area(s) with the LumaRx Full-Body Hair Removal System (female use only).

IMPORTANT: The Facial Cap has been specially designed to treat female facial hair. DO NOT use the Body Cap on the face. See section “Changing Treatment Caps” for instructions on how to exchange caps.

Where you can use the Facial Cap

The Facial Cap is intended for use on female facial hair on the cheeks, sideburns, chin, neck areas and above the lips as shown in the figure below. DO NOT use on male beards. See Contraindications on page 18.

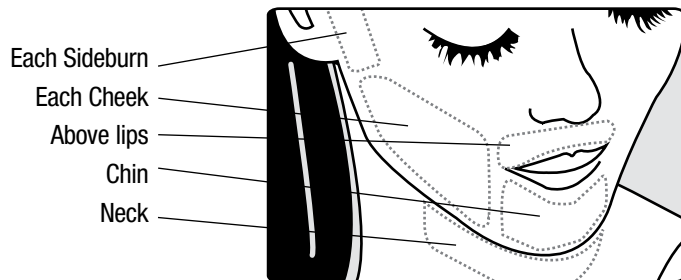

CHANGING TREATMENT CAPS

The Facial Cap is specially designed to treat Female Facial Hair, as it minimizes light leakage, and can flash every 5 seconds. The Facial Cap is marked with the icon below:



To change treatment caps, please follow the instructions provided.

1. Turn the POWER SWITCH off and unplug the POWER CORD from the electrical outlet.

Allow the system to cool down for 5 minutes before attempting to change treatment caps.

2. Hold the treatment cap with one hand while pressing both Treatment Cap Release Buttons with the other hand and pull gently. Store the unused treatment cap in a safe location (such as in the packaging) taking care not to damage the FLASH WINDOW and other components.

3. Gently push the other treatment cap back into the Hand Piece making sure both side tabs click into place. NOTE: the longer tab (with gold circuits) must be on top to properly fit.

Helpful tips for treating facial hair

- Test on a Patch of Skin.
- Use a Mirror to Guide your Treatments – using a mirror for a good view of the area to be treated and to help ensure you avoid treating the lips or any other dark spots. If the reflected light is unpleasant, you may want to close your eyes before each flash.
- Stretch Skin for Better Skin Contact – the contours of the skin may make it difficult for the skin contact sensors to be fully depressed in certain areas. You may find it easier to try and wrap your upper lip over your teeth or use your tongue to stretch your cheeks for better skin contact.

QUICK REFERENCE GUIDE

1. PREPARE SKIN (see page 27)

- Shave the hair in the treatment area
- Wipe skin clean and dry completely


2. TURN ON AND UNLOCK THE LumaRx Full-Body Hair Removal System (see page 28)




- Turn POWER SWITCH on
- Place skin in the intended treatment area on SKIN TONE SENSOR for a few seconds (TOUCHING SKIN TONE SENSOR)
- UNSUITABLE skin = “BUZZ” and orange light; the system will not unlock
- SUITABLE skin = “BEEP-BEEP,” green light, and the system will unlock
- You must unlock on each body part you intend to treat

3. SELECT ENERGY LEVEL (see page 31)

- The LumaRx Full-Body Hair Removal System starts out at the lowest level
- Test your skin at lower levels before increasing energy.
- Higher energy levels improve results but may be less comfortable
- To increase, press the ENERGY LEVEL SELECTION BUTTON
- For best results use the highest setting that is comfortable

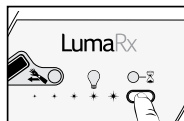

4. PLACE HAND PIECE ON AREA AND TEST A PATCH OF SKIN (see pages 28-31)

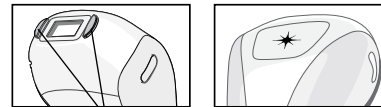

SKIN CONTACT SENSORS

- If SKIN CONTACT SENSORS are NOT in contact with skin the system will “buzz” and will not operate if you press the FLASH BUTTON. Double check that the contact sensors are in full contact with the skin.
 - If CONTACT SENSORS are in full contact with skin the system will “CHIRP” and is ready to flash.
 - If you're having trouble, try pressing the hand piece firmly against the skin.
 - See pages 28-31 for more info on testing a patch of skin.
- ### 5. TREAT YOUR SKIN (see pages 31-37)
- After you hear the “CHIRP” press the FLASH BUTTON to treat skin

- When a flash is delivered, you will:
 - See a bright flash of light
 - Hear a “POP” sound (the disposable cartridge)
 - Feel a mild sensation of warmth and tingling
 - LIFT the HAND PIECE off the skin and move it to a new treatment area
 - DO NOT flash the same treatment spot multiple times during one session
 - Try to treat in a GRID PATTERN
- Note:** DO NOT physically mark the treatment area with grid pattern as this could cause injury to the skin.
- Note:** This guide is intended to be a quick reference for experienced users of the LumaRx Full-Body Hair Removal System. If this is your first time using the LumaRx Full-Body Hair Removal System it is critical that you read the entire Use and Care Guide for complete directions for important contraindications, WARNINGS, benefits and risks and detailed use instructions.

USAGE OF THE LUMARX™ FULL-BODY HAIR REMOVAL SYSTEM (INDICATIONS)

The LumaRx Full-Body Hair Removal System is an over-the-counter device indicated for the removal of unwanted body and/or facial hair in adults with Fitzpatrick skin types I – IV.

The LumaRx Full-Body Hair Removal System is also intended for permanent reduction in unwanted hair.

Permanent hair reduction is defined as the long-term stable reduction in the number of hairs regrowing when measured at 6, 9, and 12 months after the completion of a treatment regimen.

Skin Tone

The LumaRx Full-Body Hair Removal System is intended for use only on light to medium skin tones. To determine if your Skin Tone is suitable to use with LumaRx Full-Body Hair Removal System, compare your skin in the intended treatment areas with the Skin Tone Chart. If it matches the Light to Medium Skin Tones, you can use







LumaRx Full-Body Hair Removal System. If your Skin Tone in the intended treatment area matches the Dark Skin Tones, DO NOT use the LumaRx Full-Body Hair Removal System.







The LumaRx Full-Body Hair Removal System also has a built-in SKIN TONE TESTER that is designed to help prevent flashing on skin that is too dark for treatment. See pages 28-31 for details. Test each different area you wish to treat with the Skin Tone Chart and SKIN TONE TESTER before treating. Some areas are darker than others.

If you have tanned skin in the intended treatment area, DO NOT use the LumaRx Full-Body Hair Removal System in those areas as your skin may be too dark for treatment. You should wait until your tan has completely faded from any intended treatment area before using the LumaRx Full-Body Hair Removal System.

Hair Color

The light energy targets dark coloring in hair, and therefore works best on **DARKER HAIR** (black, dark brown, brown or dark blond). **DO NOT** use the LumaRx Full-Body Hair Removal System on **LIGHTER HAIR** (light blond, red, grey or white hairs or “peach fuzz”). The LumaRx Full-Body Hair Removal System will not work on **LIGHTER HAIR**.

SKIN TONE CHART					
					
YES	YES	YES	YES	NO	NO

HAIR COLOR CHART					
Black	Brown	Dark Blond	Light Blond	Red	Peach Fuzz/ White/Gray
					
YES	YES	YES	NO	NO	NO

WHAT IS THE LUMARX™ FULL-BODY HAIR REMOVAL SYSTEM USED FOR? (INDICATIONS)


Treatment Areas

Male and Female Body Sites

The LumaRx Full-Body Hair Removal System is intended for use on male and female body hair on **underarms, bikini line, legs, arms, chest** and **stomach**. If you have unwanted DARKER HAIR in any of these areas, the illustration below lists the typical number of flashes needed to cover a 1.5" x 3" treatment area. Treat in a GRID-LIKE PATTERN – slight overlaps are OK, but DO NOT flash the same EXACT spot multiple times during one treatment session. Record treatments in the **TREATMENT CARDS** on pages 56-58.

Note: DO NOT physically mark the treatment area with grid pattern as this could cause injury to the skin.

SIZE OF FLASH WINDOW

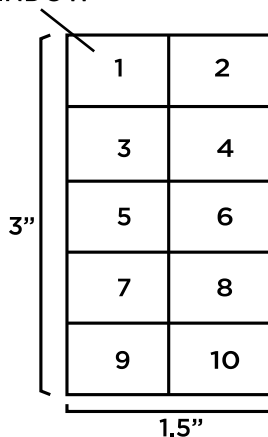

Female Facial Sites

Using the included Facial Cap, the LumaRx Full-Body Hair Removal System is intended for use on female facial hair on the **cheeks, sideburns, chin, neck areas,** and **above the lips** (see illustration). If you have unwanted DARKER HAIR in any of these areas, the illustration below lists the typical number of flashes to treat the entire site (NOTE: for spot treatments of unwanted hair, fewer flashes may be needed). Record treatments in the **TREATMENT CARDS** on pages 56-58.

DO NOT use the LumaRx Full-Body Hair Removal System on or around **EYES, EYE LASHES** or **EYE BROWS**. Doing so can cause serious and permanent eye injury.


Each Sideburn - 5 flashes

Each Cheek - 15 flashes

Above lips - 5 flashes

Chin - 10 flashes

Neck - 15 flashes

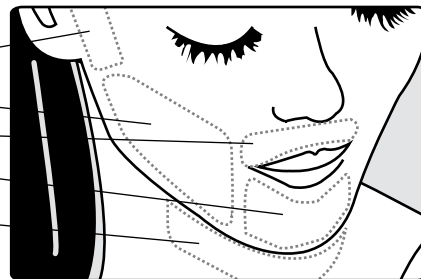

WHO SHOULDN'T USE THE LUMARX™ FULL-BODY HAIR REMOVAL SYSTEM (CONTRAINDICATIONS)

IMPORTANT SAFETY INFORMATION – READ BEFORE USE!

The LumaRx Full-Body Hair Removal System is not designed for everyone. Please read and consider the information in the following section before use.

- **DO NOT USE** the LumaRx Full-Body Hair Removal System if you have **DARK SKIN**. Darker skin may absorb too much light energy which could cause skin swelling, temporary blisters, or could change the color of your skin. Refer to the Skin Tone Chart on page 15 to determine if your Skin Tone is suitable for use with the LumaRx Full-Body Hair Removal System. You should test each different area of skin you want to treat. Some areas are darker than others.

- **DO NOT USE** the LumaRx Full-Body Hair Removal System on **TATTOOS, DARK BROWN or BLACK SPOTS** (such as dark freckles, birthmarks, moles, or warts), **LIPS, NIPPLES, GENITALS or AROUND THE ANUS**. These areas have darker skin and may absorb too much light energy which could cause skin swelling, temporary blisters, or could change the color of your skin.
- **DO NOT USE** the LumaRx Full-Body Hair Removal System **ON OR AROUND EYES, EYEBROWS or EYELASHES**. Doing so can cause serious and permanent eye injury.
- **DO NOT USE** the LumaRx Full-Body Hair Removal System anywhere on the **MALE FACE, HEAD OR NECK**.

WARNINGS

▲ The LumaRx Full-Body Hair Removal System is a powerful device. As such it should be used with special attention to safety. Please read all **WARNINGS** and safety **PRECAUTIONS** before use and follow them strictly when using the LumaRx Full-Body Hair Removal System.

▲ **DO NOT** flash in an environment where the emission of light could cause an explosion, such as in the presence of gas.

▲ **KEEP OUT OF REACH OF CHILDREN!** **DO NOT** leave the LumaRx Full-Body Hair Removal System unattended around children as the cords represent a choking hazard. Children should be supervised to ensure that they do not play with the appliance.

▲ **The LumaRx Full-Body Hair Removal System SHOULD NOT** be used by anyone under the age of 18.

▲ **DO NOT** use the LumaRx Full-Body Hair Removal System if you have a known **sensitivity to light (photosensitivity)** or are taking medications that make your skin more sensitive to light, as the LumaRx Full-Body Hair Removal System may injure your skin or make existing conditions worse.

PRECAUTIONS

The LumaRx Full-Body Hair Removal System is a powerful device. As such it should be used with special attention to safety. Please read all WARNINGS and safety PRECAUTIONS before use and follow them strictly when using the LumaRx Full-Body Hair Removal System.

- **DO NOT** use the LumaRx Full-Body Hair Removal System if you have **TANNED SKIN** in the intended treatment area, as your skin may be too dark and treatment in tanned areas could temporarily cause skin swelling, blisters, or could change the color of your skin. Wait until your tan has faded in the treatment area before treating.
- **DO NOT** use the LumaRx Full-Body Hair Removal System on skin with **MAKEUP** applied, as treating in those areas could temporarily cause skin swelling, blisters or could change the color of your skin. Make sure to thoroughly clean the intended treatment area before using the LumaRx Full-Body Hair Removal System.

- **DO NOT** flash the LumaRx Full-Body Hair Removal System on the identical treatment spot multiple times during the same session, as this could temporarily cause skin swelling, blisters, or could change the color of your skin. See instructions on pages 31-37 for more information on how to properly administer treatment.
- **DO NOT** treat any area more than once every other week, as this will not speed up results, but could increase likelihood of skin swelling, blistering or changing the color of your skin.
- **DO NOT** use the LumaRx Full-Body Hair Removal System on damaged skin, including open wounds, cold sores, sunburns, or areas where you have had a skin peel or other skin resurfacing procedure within the last 6-8 weeks.
- **DO NOT** use the LumaRx Full-Body Hair Removal System on irritated skin, including rashes or swollen skin. These conditions can make the skin more sensitive to the light treatment. Wait for

the affected area to heal before using the LumaRx Full-Body Hair Removal System.

- **DO NOT** use the LumaRx Full-Body Hair Removal System if your skin is sensitive to light. In order to determine if you can use the LumaRx Full-Body Hair Removal System, perform a test on a small patch of skin and wait 24 hours to determine if there are any adverse reactions before using the LumaRx Full-Body Hair Removal System more broadly. See steps for testing on a patch of skin on pages 27-31. Contact Customer Service (see page 51) for more information or assistance.
- **DO NOT** use the LumaRx Full-Body Hair Removal System on any area where you may someday want your hair to grow back. The results may be permanent and irreversible.
- **DO NOT** flash the LumaRx Full-Body Hair Removal System against any material other than skin.
- **DO NOT** use the LumaRx Full-Body Hair Removal System for any purpose other than hair removal.

- **DO NOT** use the LumaRx Full-Body Hair Removal System if you are or may be pregnant or are breast-feeding. The LumaRx Full-Body Hair Removal System has not been tested on pregnant women. Hormonal changes are common during pregnancy or breast-feeding and could increase skin sensitivity and the risk of skin injury.
- **DO NOT** cover any part of the LumaRx Full-Body Hair Removal System ventilation slots during use, as they help to keep the device cool during operation. Covering the ventilation slots could result in the LumaRx Full-Body Hair Removal System overheating.
- **DO NOT** use the LumaRx Full-Body Hair Removal System if the box packaging or the device appears damaged or if any part of the LumaRx Full-Body Hair Removal System body is cracked or broken. Doing so could increase the risk of skin or eye injury from the flash lamp. Contact Customer Service for assistance (see page 51).

PRECAUTIONS

- **DO NOT** allow the LumaRx Full-Body Hair Removal System to become damp or wet. **DO NOT** use the LumaRx Full-Body Hair Removal System near bathtubs, showers, basins, or other vessels containing water. Failure to follow these PRECAUTIONS can cause electrocution or electric shock. Contact Customer Service for assistance (see page 51).
- **DO NOT** operate the LumaRx Full-Body Hair Removal System with an unapproved power cord, or if it has a damaged plug or cord. Keep the LumaRx Full-Body Hair Removal System cord away from heated surfaces (such as stoves, fireplaces, radiators or furnaces). Failure to follow these PRECAUTIONS can cause electrocution or electric shock. Contact Customer Service for assistance (see page 51).
- **DO NOT** open the LumaRx Full-Body Hair Removal System (except to replace the removable TREATMENT CAP and/or LAMP CARTRIDGE) or attempt to repair your LumaRx Full-Body Hair Removal System as it contains no user-serviceable

parts. Trying to open the LumaRx Full-Body Hair Removal System may also damage the device and will void your warranty. Please contact Customer Service (see page 51) if you have a broken or damaged device. Failure to follow these PRECAUTIONS can cause electrocution or electric shock.

- Hair removal by lasers or intense pulsed light (IPL) can cause increased hair growth in some individuals.
- If you've had professional or clinical or self-administered laser or IPL or electrolysis hair removal treatments, wait at least 6 months before using the LumaRx Full-Body Hair Removal System. If you've waxed, plucked or epilated hair in the treatment area, wait at least 6 weeks before using the LumaRx Full-Body Hair Removal System. The device uses the hair under your skin to disable your hair at the root, so it is not effective if the hair has been pulled out.

RISKS

⚠ It is important to use the LumaRx Full-Body Hair Removal System properly. Make sure that you follow all the WARNINGS (see page 19), and that you do not use the device if it is not right for you. If you do not understand any of the WARNINGS, or have any questions, please contact Customer Service (see page 51).

If you use the LumaRx Full-Body Hair Removal System as directed, the risks are low:

- Your skin may become red within 24 hours after treatment. This redness generally clears up within a few hours. If the redness does not go away within 2–3 days, contact Customer Service (see page 51).
- You may develop temporary blisters, although these are unlikely for people with lighter skin tones. If you do develop a blister, stop using the LumaRx Full-Body Hair Removal System, clean the blister and cover it with a bandage and change daily to help prevent the blister from becoming infected. If you need

more information, contact Customer Service (see page 51).

- In rare cases, an outbreak of cold sores may develop in the area(s) treated that have a history of cold sores. The cold sores should heal within 5 to 7 days. If the cold sores do not heal after 7 days, contact Customer Service (see page 51).
- In rare cases, some users might see an increase in the number of hairs growing. This increase is usually temporary and due to hair growth cycles. However, if after beginning use of the LumaRx Full-Body Hair Removal System you believe that the number of hairs regrowing in the treated area has increased since you started treatment, stop using the LumaRx Full-Body Hair Removal System and contact Customer Service (see page 51).
- Skin swelling, especially around the hair follicle where you may notice the skin is slightly higher than surrounding areas. If the swelling doesn't go away within

RISKS

- 2–3 days, stop using the LumaRx Full-Body Hair Removal System and contact Customer Service (see page 51).
- Skin discoloration. If you notice that your Skin Tone in the treatment area has changed, stop using the LumaRx Full-Body Hair Removal System and contact Customer Service (see page 51).
 - Eye Injury. Flashing any laser or IPL device near or around the eyes could cause serious and permanent eye injury. This is very unlikely as long as the device is used as indicated.

WHAT TO EXPECT

During Treatment

For many people, using the LumaRx Full-Body Hair Removal System is their first experience with a light-based device. Here are some things to expect during treatment:

- A Flash of Light: The flash of light that you see is only part of the light produced by the LumaRx Full-Body Hair Removal System and should not be harmful to your eyes. The LumaRx Full-Body Hair Removal System Contact Sensors help ensure that most of the light is directed into the skin. The remaining visible light should not be harmful to the eyes, and special eye protection is not needed as long as you follow the directions for using the LumaRx Full-Body Hair Removal System in this manual. The LumaRx Full-Body Hair Removal System should NOT be used on or around the eyes, eyelashes or eyebrows.

- A Sensation of Warmth and Tingling: It is normal to feel a warm and tingling sensation with each flash. Many people describe this sensation to be like a rubber band lightly snapped against the skin. However, sensations vary by person and treatment area. If the sensation is too intense or painful, stop using the LumaRx Full-Body Hair Removal System and contact Customer Service before using again (see page 51).

Immediately After the First Treatment

- Visible changes to the hair are unlikely at this point. During and just after your LumaRx Full-Body Hair Removal System session, it is not uncommon to see some very mild redness, mostly around the hairs themselves. This is normal and typically clears up within a few hours. If the redness does not go away within 3 days, contact Customer Service (see page 51).

WHAT TO EXPECT

- If you have redness, blistering or other skin irritation in the treatment area, DO NOT perform treatment. Wait until the skin returns to normal before using the LumaRx Full-Body Hair Removal System again. If the irritation does not go away with 3 days, contact Customer Service on page 51.

One to Two Weeks after the First Treatment

- Hairs in the treated area will begin to fall out gradually, but it is unlikely you will notice significant reduction in hair at this stage.
- Make sure to complete the full treatment schedule to ensure you hit all hairs in the growth phase in the targeted area.

Two Weeks after the Last Treatment

- You should see a noticeable reduction in hairs in the treated area. Also, any new hair growth should be lighter and finer.
- At this point, additional treatments on the same area should only be done as needed to target hairs that have recently entered the active cycle. You only need to use the LumaRx Full-Body Hair Removal System for occasional touch-ups, where necessary, or to treat another area.

OPERATING INSTRUCTIONS

Once you've become familiar with the LumaRx Full-Body Hair Removal System and read all the information on the preceding pages, including the WARNINGS and PRECAUTIONS, you're ready to use the LumaRx Full-Body Hair Removal System.

This section explains step-by-step how to operate the LumaRx Full-Body Hair Removal System. Read this entire section carefully before you start treatment. Then refer back to it later as you follow the steps.

The first time you use the LumaRx Full-Body Hair Removal System, you should test it on a small patch of skin. See steps on pages 27-31.

Each time you use the LumaRx Full-Body Hair Removal System, make sure the Hand Piece, Flash Window, Skin Contact Sensor, and Skin Tone Sensor are clean and free of debris; clean the system if necessary. Refer to the User-Maintenance Instructions on pages 39-43.

STEP 1: PREPARE SKIN

Before using the LumaRx Full-Body Hair Removal System, prepare the treatment area by thoroughly Cleaning, Shaving, and Drying the intended treatment area. It is important that you remove any makeup from the treatment area.

Removing surface hair and cleaning the skin helps to minimize the sensation felt from the LumaRx Full-Body Hair Removal System flashes and allows the energy to penetrate to the hair root. Because the LumaRx Full-Body Hair Removal System works below the skin, the device is still effective even if hairs are not clearly visible.



OPERATING INSTRUCTIONS

STEP 2: TURN ON AND UNLOCK the LumaRx Full-Body Hair Removal System.

Flip the **POWER SWITCH** to turn the power on.


Each time you turn on the LumaRx Full-Body Hair Removal System, you will be prompted to test your skin tone with the SKIN TONE TESTER. The LumaRx Full-Body Hair Removal System will play a melody to remind you to test. See Figure 2.

NOTE: Using the LumaRx Full-Body Hair Removal System on darker skin tones can result in damage to the skin (such as burns, blisters or skin discoloration). To help protect against such damage, the LumaRx Full-Body Hair Removal System requires that you test your skin tone before every use. Since some skin areas are darker than others, test each different area you want to treat.


DO NOT use the LumaRx Full-Body Hair Removal System on tanned skin. Even if your natural skin tone is light, tanning makes it darker – potentially too dark for treatment. With tanned skin, wait until the tan has completely faded before using the LumaRx Full-Body Hair Removal System on that area.

Always use the Skin tone sensor to determine if your Skin is too dark to treat with the LumaRx Full-Body Hair Removal System by placing the Skin Tone Sensor on the area to be treated. Make sure your skin in this area is clean and dry, and shave it if necessary to be sure it is free of hair.

This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature.

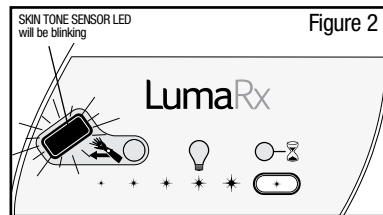


2a. Turn the **POWER SWITCH** on. The **SKIN TONE SENSOR** will flash **YELLOW**. You will hear a series of audible “BEEPS” every 5 seconds to remind you to take the test. For safety reasons, after 30 seconds of inactivity, the audible “BEEPS” will increase in frequency until the test is taken or the machine is turned off.

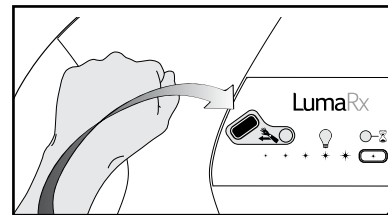

2b. Always place the skin tone sensor window directly on the area of skin you plan to treat. The reading will take place automatically.

Note: You must use the Skin Tone Sensor to test color on each body part you plan to treat.


OPERATING INSTRUCTIONS

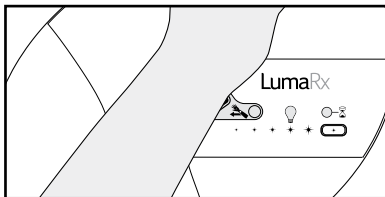


2c. Place the body part to be treated over the SKIN TONE SENSOR (touching the unit) the Sensor will immediately read the skin type.



2d. If your Skin Tone is NOT suitable for use, the SKIN TONE SENSOR turns solid YELLOW and the device will “BUZZ” once for 5 seconds. After 5 seconds, you may re-try the SKIN TONE SENSOR (step 1 above).

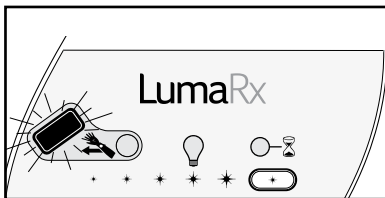

2e. If your Skin Tone is suitable for use, a high-pitched “BEEP” will sound and the SKIN TONE SENSOR will turn green and the device is ready for use.

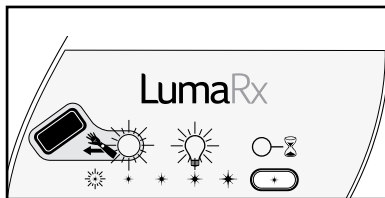

2f. Once the device is enabled, the ENERGY LEVEL LIGHT and BULB LIFE LIGHT will turn on and the internal fan will start. Your LumaRx Full-Body Hair Removal System is set to the lowest energy level setting and ready for use!

Notes:

- Some areas of skin are darker than others. You must use the SKIN TONE SENSOR on each different area you want to treat.
- If you feel the SKIN TONE SENSOR has not made a proper reading, you may try again after 5 seconds. If you believe your skin is suitable based on the Skin Tone chart on page 15 but the SKIN TONE SENSOR will not enable the device, contact Customer Service (see page 51).
- If you have tanned skin in the intended treatment area, DO NOT use the LumaRx Full-Body Hair Removal System in those areas as your skin may be too dark for treatment. You should wait until your tan has completely faded from any intended treatment area before using the LumaRx Full-Body Hair Removal System.

STEP 3: SELECT ENERGY LEVEL

Energy Level Indicator Lights determine the intensity of the LumaRx Full-Body Hair Removal System light flash delivered to your skin, from the lowest level (*) to the highest level (*****). Higher levels improve hair removal results, but may be less comfortable. To get the best results, you should use the highest setting that's comfortable.

There are 5 energy level settings on the LumaRx Full-Body Hair Removal System. Whenever the device is turned on, it is automatically set to the lowest energy level and only 1 ENERGY LEVEL INDICATOR LIGHT will be on (ENERGY LEVEL 1 in the table below). To modify the energy level, press the ENERGY LEVEL SETTING BUTTON to increase the energy level. Once the energy level has reached level 5, pressing the ENERGY LEVEL SETTING BUTTON will return the device to energy level 1. The number of ENERGY LEVEL INDICATOR LIGHTS will coincide with the change in energy level. See Figure 3 on page 32.

OPERATING INSTRUCTIONS

Figure 3

Energy Level	1	2	3	4	5
Press Energy Level Setting Button	Defaults to Level 1 upon start-up.				
Energy Level Selection Display Lights:					

3a. Upon start-up, the device sets itself to the lowest energy level (Level 1 above). The device is now ready to trigger the first flash. (See “Using the LumaRx Full-Body Hair Removal System”).

3b. To change the energy level press the ENERGY LEVEL SETTING BUTTON to increase the energy level until the desired energy level is set (for example, pressing the ENERGY LEVEL SETTING BUTTON once will increase from Level 1 to Level 2; pressing another time will increase to Level 3). ENERGY LEVEL INDICATOR LIGHTS will coincide with the energy level setting.

3c. If a higher energy level was previously selected and a lower energy level is desired, continue to press the ENERGY

LEVEL SETTING BUTTON until the desired lower energy setting is set. After reaching Level 5, pressing the ENERGY LEVEL SETTING BUTTON an additional time will reset to Energy Level 1.

As with any new skin care product, we recommend testing the LumaRx Full-Body Hair Removal System on a patch of skin before your first full treatment or before treating a new area.

STEP 4: TEST A PATCH OF SKIN

Testing the LumaRx Full-Body Hair Removal System helps you determine the appropriate Energy Level to use the first time you use the device, or for each new treatment area. Higher Energy Levels are more effective at removing hair, but may be less comfortable. To get the best results

with the LumaRx Full-Body Hair Removal System, you should use the highest Energy Level setting that’s comfortable and does not cause an adverse skin reaction.

To test the LumaRx Full-Body Hair Removal System on a new treatment area, follow these steps:

Start at Energy Level 1. The LumaRx Full-Body Hair Removal System defaults to Energy Level 1, which is where you should start the Patch of Skin Test.

Flash your Skin (see pages 34-37).

With the LumaRx Full-Body Hair Removal System at Energy Level 1, treat your skin with one single flash by firmly placing the flash window flat against skin (you should hear a “CHIRP” and the **SKIN TONE SENSOR** light should be GREEN). Press the **FLASH BUTTON** to deliver a flash.

Mild to moderate discomfort is normal. If you feel more than moderate pain, stop using the LumaRx Full-Body Hair Removal System immediately and call Customer Service (see page 51) for help.

- If a flash delivered at Energy Level 1 feels tolerable, press the ENERGY LEVEL button once to change to Energy Level 2. Move the device to a different spot and flash the LumaRx Full-Body Hair Removal System one more time. Repeat this process until you find the highest setting that feels right for you.

- Wait 24 hours and then look at the test areas. If the skin appears normal, proceed with your first full treatment at the highest setting. Record this level in your Treatment Cards on pages 56-58.

▲ If you notice blistering, burns, changes in your Skin Tone, or any other side effects, or if you are not sure about the outcome of the test, DO NOT USE the LumaRx Full-Body Hair Removal System and contact Customer Service page 51 for assistance.


Repeat the test if you decide to try the LumaRx Full-Body Hair Removal System on a different part of the face, or if your skin tone has changed (for example, from a sun tan).

OPERATING INSTRUCTIONS


STEP 5: TREAT YOUR SKIN

The LumaRx Full-Body Hair Removal System can flash approximately every 3 seconds (with the Body Treatment Cap) and is ideal for treating both larger and smaller areas.


Periodically, if either the hand piece or base unit gets too hot, the “Hourglass” icon will flash and the unit will beep. Wait for a few minutes to allow the LumaRx Full-Body Hair Removal System to cool down. There is no need to shut the unit down, as the hourglass and beeps will turn off when the LumaRx Full-Body Hair Removal System is ready to use again. The number of flashes delivered before the cool down state depends on several factors, including how fast you apply the treatments, and the room temperature. See Figure 5.

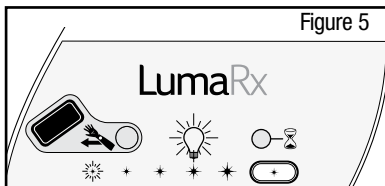



5a. Ensure that the ENERGY LEVEL is on the correct setting (See “Choose Energy Level” on page 31 for details about the different settings).

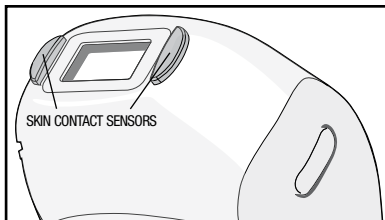

5b. Apply the treatment cap to the skin, making sure the skin is spread evenly and smoothly. Make sure the CONTACT SENSORS on the treatment cap are in full contact with the skin.

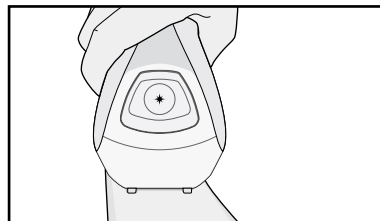


5c. If the CONTACT SENSORS are NOT in full contact with the skin, the device will not operate if you press the FLASH BUTTON and it will “BUZZ.” Double check that the contact sensors are in full contact with the skin.

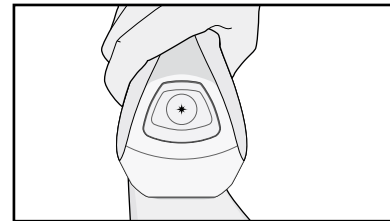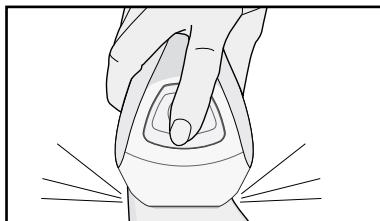

5d. If the CONTACT SENSORS on the treatment cap are in full contact with the skin, LumaRx will “CHIRP” indicating the device is OK to flash.

5e. Press the FLASH BUTTON. The device will flash a pulse of light onto your skin. You will see a flash of light and simultaneously hear a sound, which is normal for the device. You may feel a mild sensation of warmth and tingling.

OPERATING INSTRUCTIONS

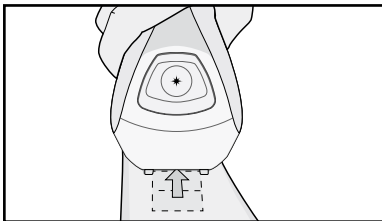

5f. Immediately after a flash is delivered, move the HAND PIECE to a new treatment location while keeping the FLASH BUTTON pressed down. After 5 seconds the BULB LIFE LIGHT will turn on again and the device will beep. If the treatment cap is in full contact with the skin and the FLASH BUTTON is pressed down, another flash will be delivered.

After a flash is delivered, lift the hand piece to the next area in a grid-format.

NOTE: If the button is kept depressed, the hand piece can be slid between spots.

About 3 seconds later (assuming Contact Sensors are in contact with skin), the device will once again emit a “BEEP” to confirm that it is ready to flash.

The device will operate if FLASH BUTTON is pressed.

5g. Try to treat in an organized, grid-like pattern to ensure you cover all hairs in the intended treatment area.
WARNING: DO NOT flash the LumaRx Full-Body Hair Removal System on the same treatment spot multiple times during the same session, as this could cause skin swelling, temporary blisters, or could change the color of your skin.

Repeat step 5 until the desired treatment area is fully covered. Depending on the area, the treatment time will vary.

NOTE: For faster coverage you may hold down the FLASH BUTTON while you lift the HAND PIECE and place it back in contact with a new area to be treated. This will cause the unit to deliver a flash to the skin as soon as full contact between the skin and the treatment cap is established.

Recommended Treatment Schedule

Your initial regimen should be spread out over 1 month to ensure that you target the most number of active hairs (including those that recently cycled from the resting phase).

We recommend that you **use the LumaRx Full-Body Hair Removal System once a week, every other week, over a month (for a total of 3 treatments)** on each area.

If you have redness, blistering or other skin irritation in the treatment area, **DO NOT** perform the treatment. Wait until the skin returns to normal before using the LumaRx Full-Body Hair Removal System again. If the irritation does not go away with 3 days, contact Customer Service on page 51.

OPERATING INSTRUCTIONS

OTHER TREATMENT TIPS

- Perform your treatments at night before going to bed, so that any lingering redness will most likely disappear by the morning.
- If you've waxed hair in the treatment area, wait 6 weeks before using the LumaRx Full-Body Hair Removal System. The device uses the hair under your skin to disable your hair at the root, so it is not effective if the hair has been pulled out.
- Shave immediately before treating the area. Doing this will let you use a higher Energy Level, because there will be less hair to heat up and the treatment will sting less. How well you shave will make a big difference in how the treatment feels. See page 27.

- Periodically, if the device gets too hot, the "Hourglass" icon will flash and the unit will beep. Wait for a few minutes to allow the LumaRx Full-Body Hair Removal System to cool down. There is no need to shut the unit down, as the hourglass icon will disappear when the LumaRx Full-Body Hair Removal System is ready to use again.
- The number of flashes delivered before the cool down state depends on several factors, including how fast you apply the treatments and the room temperature.
- Use the highest Energy Level setting that's comfortable for you and that does not cause an adverse reaction. Higher Energy Levels are more effective. See page 31.

USER-MAINTENANCE INSTRUCTIONS

When your LumaRx Full-Body Hair Removal System session is over, we suggest the following:

- Be sure to record your treatment in the Treatment Cards on pages 56-58.
- Turn the LumaRx Full-Body Hair Removal System off by pressing the POWER BUTTON.
- After each hair removal session, clean your LumaRx Full-Body Hair Removal System, especially the CAP and LIGHT WINDOW (See "Cleaning the LumaRx Full-Body Hair Removal System").
- After cleaning, store your LumaRx Full-Body Hair Removal System in its original box and keep it away from water.

CLEANING the LumaRx Full-Body Hair Removal System

After each session, clean your LumaRx Full-Body Hair Removal System, especially the LIGHT WINDOW, the CONTACT SENSORS, and the SKIN TONE SENSOR window.

⚠ DO NOT place any part of the LumaRx Full-Body Hair Removal System in water or other liquids. If any part of the LumaRx Full-Body Hair Removal System gets wet, DO NOT USE. Contact Customer Service on page 51 for help. As with any electrical device, a wet LumaRx Full-Body Hair Removal System could cause electrocution or electric shock.

- Make sure the LumaRx Full-Body Hair Removal System is turned off. Remove the POWER CORD from the device and wall socket before cleaning.
- Use a clean, lint-free cloth (one is provided with the unit) or a cotton swab to gently wipe the LumaRx Full-Body Hair Removal System surface and especially the CAP, LIGHT WINDOW, CONTACT

USER-MAINTENANCE INSTRUCTIONS

SENSORS and the SKIN TONE SENSOR window. If needed, you can dampen the cloth or cotton swab with water and mild soap. Allow the LumaRx Full-Body Hair Removal System to dry completely before using again.




- If you still need to remove dirt or debris, you can use your fingernail or a firm piece of plastic like a credit card. Be careful not to damage the FLASH WINDOW.

▲ DO NOT attempt to open or repair your LumaRx Full-Body Hair Removal System. Only authorized LumaRx Full-Body Hair Removal System repair centers are permitted to perform repairs.


REPLACING THE LAMP CARTRIDGE

The LumaRx Full-Body Hair Removal System LAMP CARTRIDGE can deliver 65,000 light flashes regardless of the energy level setting. Flash intensity is determined only according to the energy level setting of the device. There is no significant decrease of energy during the usable lifetime of the LAMP CARTRIDGE.

When a LAMP CARTRIDGE has reached 90% of its possible lifetime the BULB LIFE LIGHT on the CONTROL PANEL will turn yellow indicating that the LAMP CARTRIDGE should soon be replaced. Once all 65,000 flashes in a LAMP CARTRIDGE have been used, the BULB LIFE LIGHT will flash yellow indicating that further flashes cannot be delivered. In order to continue the hair removal session, it will now be necessary to replace the LAMP CARTRIDGE. **NOTE:** A LAMP CARTRIDGE should also be replaced if the LAMP TUBE is broken.

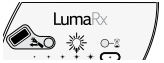
Flash Bulb Usage	0 – 90%	90 – 100%	> 100%
Flashes	0 – 58,499 flashes	58,500 – 64,999 flashes	> 65,000 flashes
Lamp Life Light			
Lamp Life Light Flash & Audible Tone	Lamp Life icon will be lit solid GREEN, with no audible tone.	Lamp Life icon will be lit solid YELLOW, with no audible tone.	Lamp Life icon will flash YELLOW, with an audible beep every 0.5 seconds until cartridge is replaced.

USER-MAINTENANCE INSTRUCTIONS

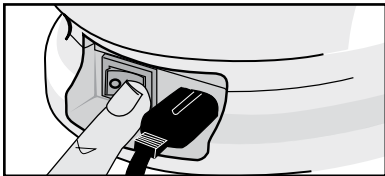


1. Turn the POWER SWITCH off and unplug the POWER CORD from the electrical outlet.

▲ Note: Allow the system to cool down for 5 minutes before attempting to change the DISPOSABLE CARTRIDGE.

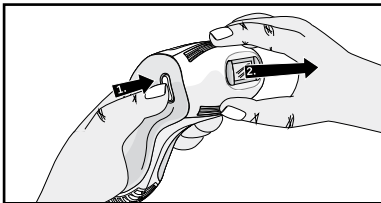


2. Grasp the treatment cap with one hand while pressing both RELEASE BUTTONS with the other hand and pull gently. Set the treatment cap in a safe location taking care not to damage the FLASH WINDOW and other components.

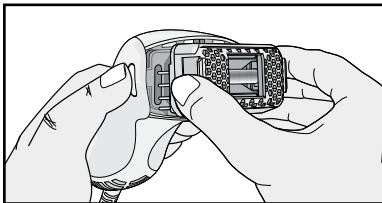


3. Grasp the LAMP CARTRIDGE and gently pull away from the HAND PIECE. Discard the used cartridge with normal household trash.

▲ WARNING: Ensure cartridge has cooled for 5 minutes before handling. Exercise care when replacing cartridge.

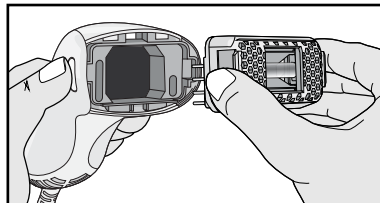

4. Push the new LAMP CARTRIDGE gently into place, taking care not to touch the LAMP TUBES or damage the pins on the cartridge. **Note:** the LAMP CARTRIDGE must be oriented so that the side with 4-pins is on the left, and the side with 1-pin is on the right.

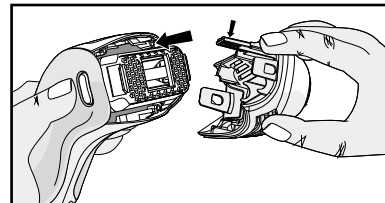

5. Gently push the treatment cap back into the HAND PIECE making sure both sides click into place. **NOTE:** the longer tab (with gold circuits) must be on top to properly fit.

DISPOSAL OF the LumaRx Full-Body Hair Removal System

Please follow the local, state and/or government requirements for disposal or recycling of electronic devices for your LumaRx Full-Body Hair Removal System, replacement head, or any of its accessories.

TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	WHAT YOU SHOULD DO
My LumaRx Full-Body Hair Removal System does not start.	The LumaRx Full-Body Hair Removal System isn't connected properly.	Make sure the power cord is properly connected to the LumaRx Full-Body Hair Removal System and that the power cord is plugged into an electrical outlet on the wall.
	The outlet isn't working.	Check if a wall switch controls the power to the wall outlet, or try plugging another device into the outlet and seeing if it works.

PROBLEM	POSSIBLE CAUSE	WHAT YOU SHOULD DO
When I place the LumaRx Full-Body Hair Removal System against my skin, nothing happens. The HAND PIECE doesn't beep, even after a few seconds.	The LumaRx Full-Body Hair Removal System is locked.	Unlock the LumaRx Full-Body Hair Removal System with the SKIN TONE SENSOR, based on steps on pages 28-31. Note that the LumaRx Full-Body Hair Removal System must be unlocked again if you don't use it for 5 minutes. If you take a break in the middle of a treatment, or turn off the device, you will need to unlock the LumaRx Full-Body Hair Removal System again to use it. See steps on pages 28-31 for "Test a Patch of Skin" using the Skin Tone Sensor.
	The CONTACT SENSORS aren't in full contact with the skin.	Check that both CONTACT SENSORS are touching your skin. The LumaRx Full-Body Hair Removal System won't work unless it's flat against your skin. See pages 34-35 for this situation.
	You are trying to treat a small area, such as fingers or toes.	You can try pulling your skin to make a bigger flat area, pushing part of your skin so that the window can reach it (such as another finger) next to the part you want to treat. If you can't get the entire window to touch your skin, the LumaRx Full-Body Hair Removal System will not work.


TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	WHAT YOU SHOULD DO
A light flash is not emitted when I press the FLASH BUTTON.	The CONTACT SENSORS are not in full contact with the skin.	Check that both CONTACT SENSORS are touching your skin. The LumaRx Full-Body Hair Removal System won't work unless it's flat against your skin. See pages 34-35 for this situation.
	You have not waited enough time between flashes for the LumaRx Full-Body Hair Removal System bulb to recover.	The LumaRx Full-Body Hair Removal System requires about 3-4 seconds to recover between flashes. When your LumaRx Full-Body Hair Removal System is ready to flash (and assuming the CONTACT SENSORS are flush against the skin), the device will "CHIRP" to confirm that it is ready to flash. See pages 34-37 for this situation.
	The LAMP CARTRIDGE has flashed 65,000 times and needs to be replaced. You should see the Light Bulb flashing YELLOW.	Make sure that the BULB LIFE LIGHT is either solid green or solid yellow. If it is flashing yellow, disconnect the LumaRx Full-Body Hair Removal System and replace the LAMP CARTRIDGE. Refer to "Replacing the LAMP CARTRIDGE" on page 41.

PROBLEM	POSSIBLE CAUSE	WHAT YOU SHOULD DO
The "hourglass" icon is flashing yellow and the LumaRx Full-Body Hair Removal System is beeping	The LumaRx Full-Body Hair Removal System hand piece needs some time to cool down.	Wait for a few minutes for the LumaRx Full-Body Hair Removal System to cool down. There is no need to shut the unit down, as the hourglass icon and beeps will turn off when the LumaRx Full-Body Hair Removal System is ready to use again.

TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	WHAT YOU SHOULD DO
<p>The outer ENERGY LEVEL lights are flashing alternately and the LumaRx Full-Body Hair Removal System is beeping every ½ second.</p> <p>⊛ * * * *</p>	<p>The LumaRx Full-Body Hair Removal System isn't working properly.</p>	<p>Contact Customer Service (see page 51) for help.</p>

PROBLEM	POSSIBLE CAUSE	WHAT YOU SHOULD DO
<p>All the lights on the system are flashing and the unit is beeping constantly.</p> 	<p>The treatment cap is removed or not properly attached.</p>	<p>Check to ensure that the treatment cap is properly attached to the HAND PIECE. See “REPLACING THE LAMP CARTRIDGE” on page 41 for instructions on how to check that the treatment cap is properly attached. When the treatment cap is properly attached, you should hear a “click” and the lights and audible tones will stop.</p>
<p>I used the SKIN TONE SENSOR and it showed an ORANGE light. I can't unlock the system.</p>	<p>Your skin is too dark (see the Skin Tone Chart on page 15).</p>	<p>You cannot use the system with your Skin Tone. Please contact Customer Service (see page 51) for help.</p>
	<p>You didn't use the SKIN TONE SENSOR correctly.</p>	<p>Try using the SKIN TONE SENSOR again – refer to pages 28-31 for instructions on using the SKIN TONE SENSOR.</p>
	<p>The SKIN TONE SENSOR isn't working.</p>	<p>Please contact Customer Service (see page 51) for help.</p>

TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	WHAT YOU SHOULD DO
Dirt or debris is stuck to the window.	The LumaRx Full-Body Hair Removal System needs to be cleaned.	See pages 39-40 for instructions on how to clean your LumaRx Full-Body Hair Removal System. If you still can't get it clean, contact Customer Service on page 51 for help.

▲ DO NOT attempt to open or repair your LumaRx Full-Body Hair Removal System. Only authorized LumaRx repair centers are permitted to perform repairs.

Opening the LumaRx Full-Body Hair Removal System may expose you to dangerous electrical components and to pulsed light energy, either of which may

cause serious bodily damage and/or permanent eye injury. Trying to open the LumaRx Full-Body Hair Removal System may also damage the device and will void your warranty. Please contact LumaRx Customer Service on page 51 if you have a broken or damaged device in need of repair.

CONTACT INFORMATION

Customer Service: 800-736-4648

Online: www.Luma-Rx.com

Address: Spectrum Brands, Inc.
3001 Deming Way
Middleton, WI 53562

SPECIFICATIONS

Spot Size	3.0cm ² (for Body Cap), 2.0cm ² (for Facial Cap)
Speed	3 seconds (for Body Cap), 4 seconds (for Facial Cap)
Technology	Xenon Flash Lamp (IPL)
Max Energy Level	8 J/cm ² ($\pm 20\%$)
Wavelength	650-1200 nm (for Body Cap), 600-1200nm (for Facial Cap)
Electrical Requirements	120-230VAC, 1.8A
Operation and Safety	UL 60335, EN60335
Package Size	10.125" x 9.625" x 8.3125"

CUSTOMER SERVICE

For more information about LumaRx, please visit www.Luma-Rx.com. If your LumaRx Full-Body Hair Removal System is broken, damaged, in need of repair, or for any other LumaRx Full-Body Hair Removal System user assistance, please contact LumaRx Customer Service: 1-800-736-4648.

PERFORMANCE GUARANTEE

Limited Two-Year Warranty

Spectrum Brands, Inc. warrants this product against any defects that are due to faulty material or workmanship for a two-year period from the original date of consumer purchase. This warranty does not include damage to the product resulting from accident or misuse.

If the product should become defective within the warranty period, we will replace it free of charge. Return your product and sales receipt with your name, address and day time phone number to: Spectrum Brands, Returns Center, 507 Stokely Dr., P.O. Box 1, Deforest, WI 53532. For more information call 800-736-4648.

KEEP ORIGINAL SALES RECEIPT AS PROOF OF PURCHASE FOR WARRANTY PURPOSES.

This warranty does not cover products damaged by the following:

- Accident, misuse, abuse or alteration of the product
- Servicing by unauthorized persons
- Use with unauthorized accessories
- Connecting it to incorrect current and voltage
- Any other conditions beyond our control

PERFORMANCE GUARANTEE

SPECTRUM BRANDS, INC. SHALL NOT BE RESPONSIBLE FOR ANY INCIDENTAL, SPECIAL OR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT. ALL IMPLIED WARRANTIES, INCLUDING BUT NOT LIMITED TO IMPLIED WARRANTIES OF FITNESS AND MERCHANTABILITY, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE.

This warranty gives you specific legal rights and you may also have other rights which vary from state to state. Some states do not allow the exclusion or limitation of incidental, special or consequential damages.

Performance Guarantee

Spectrum Brands, Inc. Money-Back Offer
If within 90 days after you have purchased your LumaRx branded product, you are not satisfied and would like a refund, return it with the sales slip indicating purchase price and date of purchase to the retailer from whom it was purchased. Spectrum Brands, Inc. will reimburse all retailers who accept the product within 90 days from the date of purchase.

If you have any questions concerning the money-back guarantee, please call 800-736-4648 in the U.S.

120-230 VAC, 50/60 Hz

® Registered Trademark of Spectrum Brands, Inc., (SBI)

Questions or comments:

For more information please visit www.Luma-Rx.com. If your LumaRx Full-Body Hair Removal System is broken, damaged, in need of repair, or for any other user assistance, please contact Customer Service: 800-736-4648 or visit www.Luma-Rx.com

Dist. by:
Spectrum Brands, Inc.,
Middleton, WI 53562
01/15

Made in China.

T22-0003041
610-00019, Rev D

© 2014 SBI



TREATMENT CARDS

For best results, you need to adhere to the prescribed care regimen. Use the **Suggested Number of Flashes** below for each treatment site, to ensure proper treatment. Fill out the Treatment Cards below to help keep track of your treatments.


Body Hair

For a typical 1.5" x 3" treatment area, you should use 10 flashes to fully treat that complete area.

Treat in an organized grid pattern – slight overlaps are OK, but do NOT treat the same spot multiple times during one session.

Note: DO NOT physically mark the treatment area with grid pattern as this could cause injury to the skin.

SIZE OF
FLASH
WINDOW

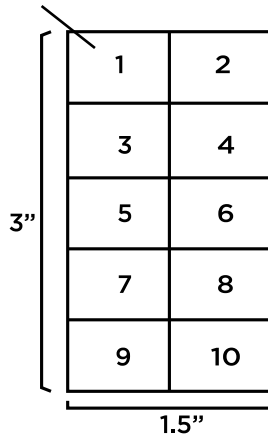

Recommended Treatment Schedule:

- Use LumaRx Full-Body Hair Removal System once every other week (total of 3 treatments) on each area.
- If you have redness, blistering or other skin irritation in the treatment area, DO NOT perform treatment. Wait until the skin returns to normal before using LumaRx Full-Body Hair Removal System again. If the irritation does not go away within 3 days, contact Customer Service on page 51.

TREATMENT CARDS

Treatment Area: Above the Lips				
Typical User: 5 flashes over entire area				
#	Date*	Areas Treated	Energy Level	# of Flashes
1				
2				
3				
4				
5				
6				
7				
8				

Treatment Area: Above the Lips				
Typical User: 5 flashes over entire area				
#	Date*	Areas Treated	Energy Level	# of Flashes
1				
2				
3				
4				
5				
6				
7				
8				

Treatment Area: Each Cheek				
Typical User: 15 flashes over entire area				
#	Date*	Areas Treated	Energy Level	# of Flashes
1				
2				
3				
4				
5				
6				
7				
8				

Treatment Area: Each Cheek				
Typical User: 15 flashes over entire area				
#	Date*	Areas Treated	Energy Level	# of Flashes
1				
2				
3				
4				
5				
6				
7				
8				

Treatment Area: Each Sideburn				
Typical User: 5 flashes over entire area				
#	Date*	Areas Treated	Energy Level	# of Flashes
1				
2				
3				
4				
5				
6				
7				
8				

Treatment Area: Each Sideburn				
Typical User: 5 flashes over entire area				
#	Date*	Areas Treated	Energy Level	# of Flashes
1				
2				
3				
4				
5				
6				
7				
8				

Treatment Area: Chin				
Typical User: 10 flashes over entire area				
#	Date*	Areas Treated	Energy Level	# of Flashes
1				
2				
3				
4				
5				
6				
7				
8				

Treatment Area: Chin				
Typical User: 10 flashes over entire area				
#	Date*	Areas Treated	Energy Level	# of Flashes
1				
2				
3				
4				
5				
6				
7				
8				

TREATMENT CARDS

Treatment Area: Neck				
Typical User: 15 flashes over entire area				
#	Date*	Areas Treated	Energy Level	# of Flashes
1				
2				
3				
4				
5				
6				
7				
8				

Treatment Area: Neck				
Typical User: 15 flashes over entire area				
#	Date*	Areas Treated	Energy Level	# of Flashes
1				
2				
3				
4				
5				
6				
7				
8				

Treatment Area:				
Typical User: 10 flashes over 1.5" x 3" area				
#	Date*	Areas Treated	Energy Level	# of Flashes
1				
2				
3				
4				
5				
6				
7				
8				

Treatment Area:				
Typical User: 10 flashes over 1.5" x 3" area				
#	Date*	Areas Treated	Energy Level	# of Flashes
1				
2				
3				
4				
5				
6				
7				
8				