
powered by **BRAUN**

PRO 1500-4000  
Deep Sweep™ 1500-4000

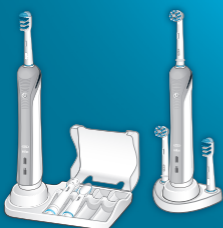English 3

Español 10

Français 17











**Charger Type 3757**  
**Cargador modelo 3757**  
**Type de chargeur 3757**

**Handle Type 3756**  
**Mango modelo 3756**  
**Type de manche 3756**

**Internet:**  
[www.oralb.com](http://www.oralb.com)  
[www.oralb.ca](http://www.oralb.ca)  
[www.braun.com](http://www.braun.com)  
[www.service.braun.com](http://www.service.braun.com)

97586353/IX-15  
EN/E/F

Braun GmbH  
Frankfurter Strasse 145  
61476 Kronberg  
Germany/Alemania/Allemagne

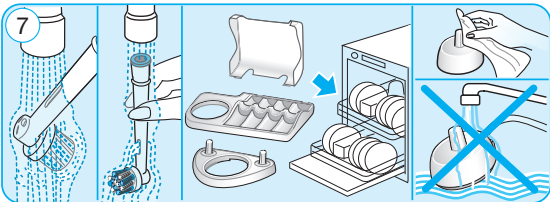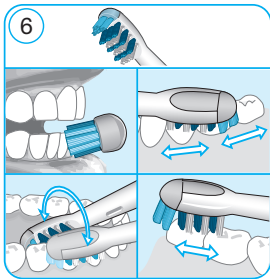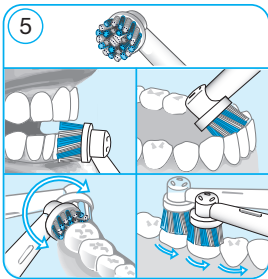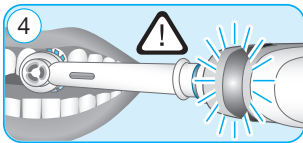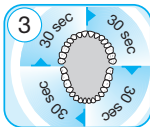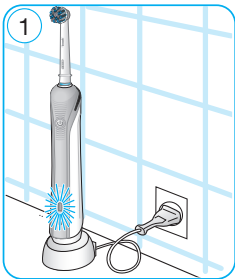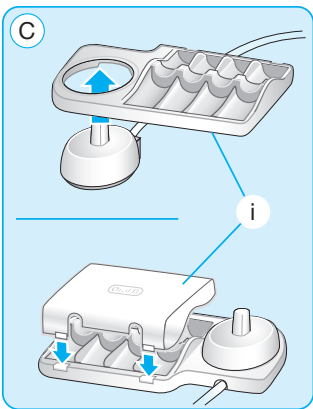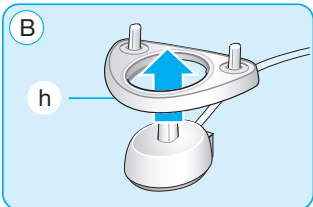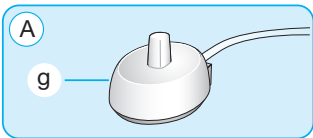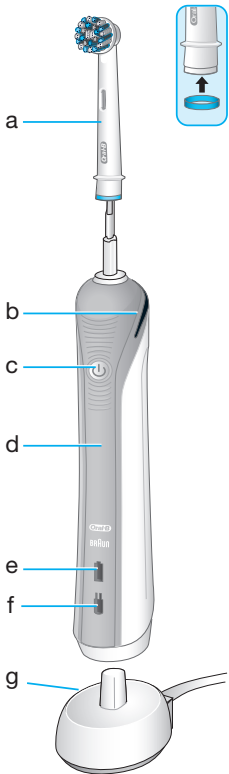---

# IMPORTANT SAFEGUARDS

Your **Oral-B®** toothbrush has been carefully designed to offer you and your family a unique brushing experience that is both safe and effective.

Read all instructions before using.

## DANGER

In order to reduce the risk of electrocution:

1. Do not place or store the charger where it can fall or be pulled into a tub or sink.
2. Do not place the charger in water or other liquid.
3. Do not reach for a charger that has fallen into water. Unplug immediately.
4. Do not use while bathing.

## WARNING

To reduce the risk of burns, electrocution, fire or injury:

1. Close supervision is necessary when this product is used by, on, or near children or invalids.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Keep cord away from heated surfaces.
4. Never drop the appliance or insert any object into any opening of the appliance.
5. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged or if it has been dropped into water. If any of these occur, please contact a service center and return the product for examination and repair.

## IMPORTANT

1. Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral-B Service Center. A damaged or non-functioning unit should no longer be used.
2. If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible.
3. Do not touch the power plug with wet hands. This can cause an electric shock.
4. When unplugging, always hold the power plug instead of the cord.
5. Do not modify or repair the product. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Service Center. Do not disassemble the product except when disposing of the battery.
6. When taking out the battery for disposal, use caution not to short the positive (+) and negative (-) terminals.
7. This product is not intended for use by children under age three.
8. This toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
9. Replace brush (head) every 3 months or sooner if brush head becomes worn.

## CAUTION

If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

## SAVE THESE INSTRUCTIONS

---

# English

## Specifications

For voltage specifications refer to the bottom of the charging unit.

## Description

- a Brush head
- b Pressure sensor light (depending on model)
- c On/off button (mode selection)
- d Handle
- e Charge indicator light
- f Low charge indicator light
- g Basic charger
- h Brush head holder (depending on model)
- i Brush head compartment with protective cover (depending on model)

## Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the basic charger (g) or the assembled charging unit (picture B/C, page 2) into an electrical outlet and place the toothbrush handle (d) on the charging unit.
- The green charge indicator light (e) flashes while the handle is being charged. Once it is fully charged, the light turns off. A full charge can take up to 22 hours and enables up to seven days of regular brushing (twice a day, 2 minutes) (picture 1, page 2).  
Note: After deep discharge, the green light might not flash immediately; it can take up to 10–15 minutes.
- When the rechargeable battery is low, the red low charge indicator light (f) starts flashing after turning off the toothbrush.
- For everyday use, the toothbrush handle can be stored on the charging unit to maintain it at full power. Overcharging is impossible.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

## Using your toothbrush

Wet brush head and apply any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching on the appliance (picture 2, page 2). When brushing your teeth with one of the **Oral-B oscillating-rotating brush heads** guide the brush head slowly from tooth to tooth, spending a few seconds on each tooth surface (picture 5, page 2). When using the Oral-B **«Deep Sweep» brush head** place the toothbrush bristles against the teeth at a slight angle towards the gumline. Apply light pressure and start brushing in back and forth motions, just like you would do with a manual toothbrush (picture 6, page 2).

With any brush head start brushing the outsides, then the insides and finally the chewing surfaces. Brush all four quadrants of your mouth equally. You may also consult your dentist or dental hygienist about the right technique for you.

In the first few days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist. If you have sensitive teeth and gums, Oral-B recommends using the «Sensitive» mode (optional in combination with the Oral-B Sensitive Gum Care brush head).

## Brushing modes (depending on model)

1-Mode	2-Mode	3-Mode	Brushing modes (for handle type 3756)
✓	✓	✓	<b>Daily Clean</b> – Standard for daily mouth cleaning
	✓	✓	<b>Sensitive</b> – Gentle, yet thorough cleaning for sensitive areas
		✓	<b>Whitening</b> – Exceptional polishing for occasional or everyday use

How to switch between modes:

Your toothbrush automatically starts in «Daily Clean» mode.

To switch to other modes, successively press the on/off button.

If you wish to turn off your toothbrush from «Sensitive» mode e.g., push and hold down the on/off button until the motor stops.

### Pressure Sensor

Helping protect your teeth and gums from hard brushing, your toothbrush has a pressure control feature installed. If too much pressure is applied, the red pressure sensor light (b) will light reminding you to reduce pressure (picture 4, page 2) (depending on model - starting with model 2000).

In addition the movement of the brush head will continue but its pulsation will stop and you will also hear a different sound while brushing. Periodically check the operation of the pressure sensor by pressing lightly on the brush head during use.

### Professional Timer

A short stuttering sound at 30-second intervals reminds you to brush all four quadrants of your mouth equally (picture 3, page 2). A long stuttering sound indicates the end of the professionally recommended 2-minute brushing time. The elapsed brushing time is memorized when the handle is briefly switched off during brushing. When pausing longer than 30 seconds the timer resets.

## Brush heads

Oral-B offers you a variety of different Oral-B brush heads that fit your Oral-B toothbrush handle.

Our **oscillating-rotating brush heads** can be used for tooth-by-tooth precision cleaning.


#### Oral-B CrossAction brush head

Our most advanced brush head. Angled bristles for a precise clean. Lifts and powers away plaque.


#### Oral-B Precision Clean brush head

provides effective and gentle cleaning action for both teeth and gums. The brush head cups each tooth with longer and denser bristles and reach deep between teeth at the same time.
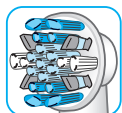
#### Oral-B FlossAction brush head

features micropulse bristles that allow for an outstanding plaque removal in the interdental areas.
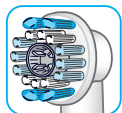
#### Oral-B 3D White brush head

features a special polishing cup to naturally whiten your teeth. Please note that children under 12 years old should not use the Oral-B 3D White brush head.


#### Oral-B Sensitive Gum Care brush head

provides thorough cleaning for sensitive teeth and gums. It has an extra soft bristle texture that is gentle on teeth and gums.

#### Oral-B Ortho brush head


with special tufting pattern thoroughly cleans and gently removes plaque from teeth with braces.
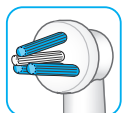
#### Oral-B Power Tip brush head

is designed for cleaning between teeth as well as around bridges, crowns and implants. Position the Power Tip at a 45 degree angle to the gumline. Slowly guide the Power Tip along the gumline, pausing between teeth for excellent results.

The **dynamic «Deep Sweep»** brush head can be used with all common brushing techniques.

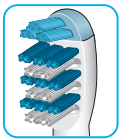Oral-B Deep Sweep brush head features triple zone cleaning action for outstanding plaque removal, even in between teeth.

Oral-B brush heads feature light blue INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the blue color will fade halfway down approximately in a 3-month period. If the bristles splay before the color recedes, you may exert too much pressure on teeth and gums.

We do not recommend using the Oral-B FlossAction or Oral-B 3D White brush head with braces. You may use the Oral-B Ortho brush head, specifically designed to clean around braces and wires.

## Cleaning recommendations

- After brushing, rinse your brush head carefully under running water for some seconds with the handle switched on.
- Switch off the handle and remove the brush head.
- Clean both parts separately under running water then wipe them dry before reassembling the toothbrush.
- Disassemble charging unit before cleaning. Never place the charger in dishwasher or water; it should be cleaned with a damp cloth only (picture 7, page 2).
- Brush head holder (h), brush head compartment and protective cover (i) are dishwasher-safe.

## Environmental Notice


Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

## Removing batteries

### Handle battery

At the end of the product's useful life, open the handle as shown, remove the rechargeable battery and dispose of it according to local environmental regulations.

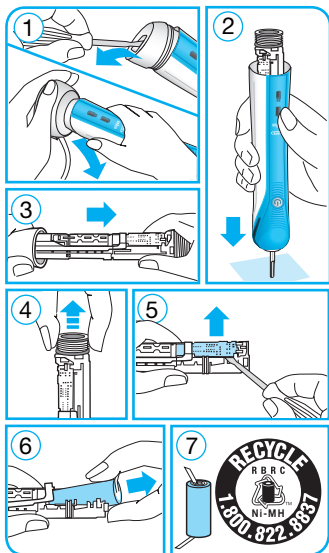**Caution! Opening the handle will destroy the appliance and invalidate the warranty.**

Subject to change without notice.

## Satisfaction Guaranteed or Your Money Back!

Try Oral-B for 60 days, starting from the day of purchase. If you are not 100% satisfied, receive a full refund via pre-paid card. No questions asked. Please return your product along with your name, address, contact phone number and the original cash register receipt indicating purchase price and date to the address below to receive a pre-paid card in the amount of your purchase (refunds in Canada include tax and Canada Post 1st class mail postage).

### **U.S. only:**

Oral-B 60 Day Satisfaction Guarantee  
4608 Appliance Drive  
Suite 2  
Belcamp, MD 21017

### **Canada only:**

Attention: 60-Day Challenge  
Braun MBG  
PO Box 12060  
Saint John, NB  
E2L 5E7

Your Oral-B toothbrush must be returned in its original packaging and be postmarked no later than 60 days after the date of purchase. Please allow 6–8 weeks for handling. In the U.S., the shipment cost to return the product is the responsibility of the consumer.

Sorry: No refunds will be made for any Oral-B product damaged by accident, neglect, unreasonable use or lost in the post.

Terms and conditions apply. For further details please contact our customer service

### **Warranty replacement brush heads**

Oral-B warranty will be voided if electric rechargeable handle damage is found to be attributed to the use of non-Oral-B replacement brush heads.

Oral-B does not recommend the use of non-Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

---

# FOR USA ONLY

## **Limited 2-Year Warranty**

Oral-B warrants this appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase at retail.

If the appliance exhibits such a defect, Oral-B will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an authorized Oral-B Service Center. Proof of date of original purchase is required.

The cost to ship the product to an authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Oral-B. This warranty does not cover finishes or normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL ORAL-B BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you. This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state.

## **How to obtain service in the USA**

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

*To obtain service:*

- A. Carry the product in to the authorized Oral-B Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

*To obtain original Oral-B replacement parts and accessories:*  
Contact an authorized Oral-B Service Center.  
Please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

*To order original Oral-B replacement parts and accessories (including replacement brush heads) by phone:*  
Please call 1-800-566-7252.

*To speak to an Oral-B Consumer Service Representative:*  
Please call 1-800-566-7252.

---

# FOR CANADA ONLY

## **Limited 2-Year Warranty**

In the event this appliance fails to function within the specified warranty period because of defects in material or workmanship, and the consumer returns the unit to an authorized service center, Braun Canada will, at its option either repair or replace the unit without additional charge to the consumer. This warranty does not cover any product which has been damaged by dropping, tampering, servicing performed or attempted by unauthorized service agencies, misuse or abuse. This warranty excludes units which have been modified or used for commercial purposes.

To the extent allowed by law, Braun Canada shall not be responsible for loss of the product, loss of time, inconvenience, commercial loss, special or consequential damages. The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

## **How to obtain service in Canada**

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

### *To obtain service:*

- A. Carry the product in to the authorized Oral-B Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

### *To obtain original Oral-B replacement parts & accessories:*

Contact an authorized Oral-B Service Center.

Please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

### *To speak to an Oral-B Service Representative:*

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