



For more information or replacement brushes  
visit us at [taoclean.com](http://taoclean.com)



Call or text: +1-218-382-5326  
Southern California, USA



Email: [info@taoclean.com](mailto:info@taoclean.com)  
Southern California, USA


# AURA CLEAN SYSTEM

## Sonic Toothbrush & Cleaning Station

### USER MANUAL



BA-0101

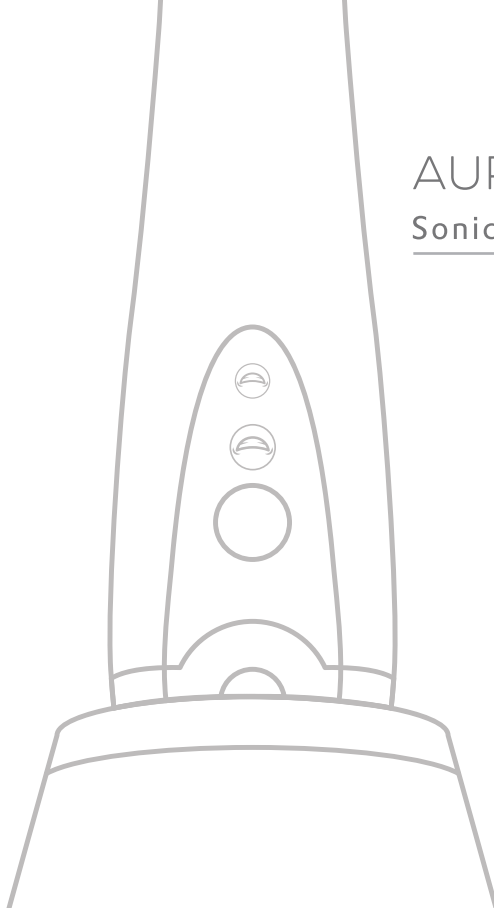

# Contents

Introduction	03	Dashboard	16
Some very important notes before use	04	Activating the brush	17
What's in the box?	08	Brushing modes	17
Charging mode	09	Brushing technique	18
Cleaning mode	10	Brush maintenance	20
Drying mode	11	Base station cleaning	21
Brush assembly	12	Frequently Asked Questions	22
Connecting power	13	Warranty	25
Brush docking	14	Contact	27
Low battery	15		

## AURA CLEAN SYSTEM WELCOME TO A NEW ERA OF CLEANLINESS


If you are reading this fresh from the box of your Aura Clean System, we salute you as a smart cookie. We would like to get you up and running as fast as possible. To ensure the best from your Aura Clean System, please take the short time needed to read these instructions.


The Aura Clean System has been lovingly designed and made to serve you loyally everyday to help you keep to your meticulous cleanliness standards.


Treat reading this manual like our first date. We'll establish some important ground rules, then we'll whisk you through some very easy to follow instructions to get the best out of our intimate relationship. We already acknowledged that you must be a smart cookie, so please excuse the formality of introducing some very necessary safeguards before we get underway.


# Attention Please!

## When using electrical devices, basic safety precautions should always be followed

 Be safe and **DO NOT** reach for an Aura Clean Base Station that has fallen into water. First unplug immediately, then gather the Aura Clean Base Station and wait until fully dry before activating.

 Never force the Aura Clean Base Station plug into an outlet. If it does not easily fit into the outlet, discontinue use and return to your vendor or TAO Clean for service.

 Please **DO NOT** plug the Aura Clean Base Station into an outlet with a voltage or frequency other than specified on the charger. Please be aware voltage converters **DO NOT** guarantee voltage compatibility. So be vigilant and most of all ...be safe.

 Now this seems obvious, never use the Aura Clean Base Station charger if it has a damaged cord or plug adaptor. Please keep the cord away from heated surfaces and do not use the power adaptor if dropped into water.



The Aura Clean System can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they understand the hazards involved and have been given supervision or instruction concerning use of the Aura Clean System in a safe way. Phew.



Please ensure that children do not play with the Aura Clean System. It is not a toy or real rocket and space station. It's for cleaning intensively and sensibly. This goes for adults too.



In case of any mechanical failure within the warranty period, please return your Aura Clean System for service. Do not attempt to dismantle.



### WARNING

This appliance contains a UV emitter. Do not stare at the light source.

# Important!

## The Aura Clean Brushing Experience

When our design council first gathered within the marbled halls of TAO Clean, situated high atop Mt. Cleanliness, it was decreed that every object we create must deliver a wonderful user experience from start to finish. It was also decreed that all our life tools must be built on outstanding aesthetic design and be both easy to adopt and seamlessly entwined with your daily routine. These words have been our mantra to deliver you the Aura Clean System Experience. Now, in order to get the best of this intention, please be aware of the following:

This Aura Clean brush is designed to clean your teeth, gums and tongue only. Use this device only for its intended use. This may seem super obvious, however some people have an insatiable curiosity.

Please discontinue use of the Aura Clean System and contact a physician/dentist if you experience discomfort or pain. The Aura Clean System is your partner for helping you achieve optimum oral wellness in the cleanest way possible. In fact your Aura Clean System took a vow never to willfully hurt you before it was carefully nestled into its delivery crib. It would never want to hurt you. Ever. Never ever ...ever.

Now, this is also very important. Please consult your dental professional before you use this device if you have had oral or gum surgery in the previous 2 months.

Contact your dental professional if excessive bleeding occurs after using the Aura Clean Brush or bleeding continues to occur after 1 week of use. You may be pressing too hard so please look at page 18 for brushing tips.

Your Aura Clean System has been tested and is compliant with safety standards for electromagnetic devices.

In any case, please consult your physician prior to using the Aura Clean System if you have medical concerns. We want to keep you clean. We want to keep you healthy. Most importantly we want to keep you that way.

Do not use attachments such as power cords or brush heads other than those recommended by TAO Clean. Otherwise your Aura Clean System will become very sad. In other words, it won't work as we've designed it to work.

If this Aura Clean Brush is dropped, don't worry. It won't blame you. However, the brush head should be immediately replaced before the next use; even if there is no visible damage to the brush head.

# What's in the Box?

Carefully open the delivery crib for your Aura Clean System. There it is - nestled inside asleep and ready for awakening. It's probably dreaming about you at this very moment and looking forward to meeting its partner. Let's check and make sure you have everything.

Aura Clean System check list:

- 1 Power Adaptor & Cable
- 2 Aura Clean Base Station
- 3 Aura Clean Brushes
- 4 Aura Clean Brush Handle



# Charging Mode

So you can take it for a spin right away, your Aura Clean Brush Handle should have enough charge to use at least once. Afterwards you will need to charge the battery.

**Note:** The battery needs 24 hours to fully charge the first time.

This could never be simpler (unless you've ignored connecting the power to the Aura Clean Base Station).

**Just place the Brush Handle in the plugged-in Base Station and it will automatically start to charge.**

Placing your Brush Handle in the plugged-in Base Station will immediately start the cleaning cycle (see page 10) and at the same time start the charging process. After the Cleaning process has finished (approximately 5 minutes) the

badge light will change from blue to white. It is now feeding the Aura Clean Brush with all that clean electricity. Yum.

Keeping the Brush Handle on its Base Station not only ensures the brush head is protected in a nice clean environment, it will also keep it fully charged and ready for use.

If possible, try to keep your Aura Clean Brush Handle on the base when not in use.



# Cleaning Mode


Our clever design of the Aura Clean Base Station protects and encloses the Aura Clean Brush from environmental pollutants.

Cleaning the Brush takes place within the Base Station core. Here lives the magical Aura Clean Halo. The Aura Clean Halo is a circular UV lamp central to the cleaning and sanitizing cycle. The intense 360 degree light source gives germs nowhere to hide. If you must know, the lamp's UV-C rays destroy the DNA of pathogens keeping your brush from turning into some ugly biology experiment.

Just drop the Brush in the Base Station and let it do its thing.

As the cleaning cycle begins the LED light will change from white to blue and slowly pulse like a mini germ-killing reactor...

After the cycle is finished the LED light will change back to white.

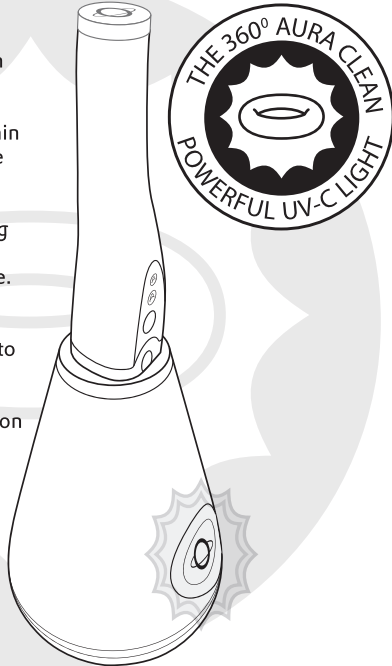

# Drying Mode

**Killing 99.9% of germs is futile if the brush is kept wet.**


That 0.1% that remains has a nasty habit of making a comeback.

To combat this, the Aura Clean Base Station has a dryer that drives water from the brush.


In fact the Base Station dries the toothbrush up to six times faster than regular UV sanitizer chambers.



# Brush Assembly

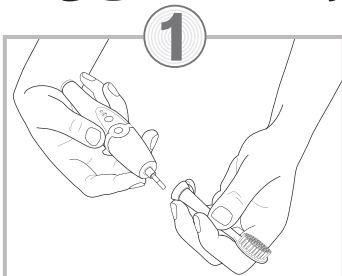

1  
Connect the brush to the brush handle. Locate hole on stem and ...

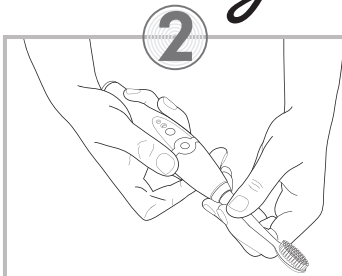


2  
Gently and firmly lock into place making sure the brush bristles are facing the same direction as the controls on the brush handle.

## Brush Replacement

Replace Aura Clean brush heads every 3 months to achieve the best results or change before if you feel you need to. Some people like it that way.


Use only official Aura Clean replacement Brush Heads. Fakes won't get the job done.

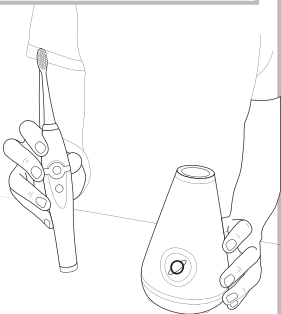

# Connecting Power


Time to feed your Aura Clean Base Station and wake it up to perform functions. Carefully plug the power adaptor into the small socket located below the waistline at the back of the Base Station.

Use the correct orientation of the plug into the socket. Now you are ready to plug the adaptor into an electrical outlet.

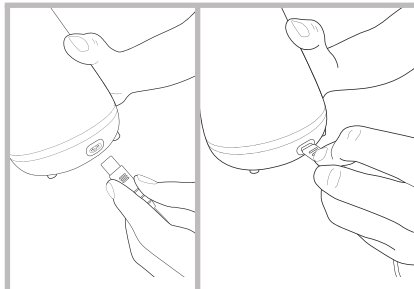

# Attention!



To ensure your Aura Clean Base Station performs its tasks perfectly, please make sure your Base Station is placed on a clean, flat horizontal surface.

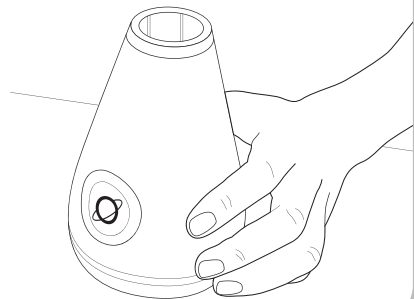

# Brush Docking

Gently place your brush into the Base Station. The Brush fits to the Base Station perfectly. Take care in locating the Brush so that it drops into the Base effortlessly.

**Hint:** the control panel on the brush should line up with the badge on the Base. Did we mention to do this carefully? Good.

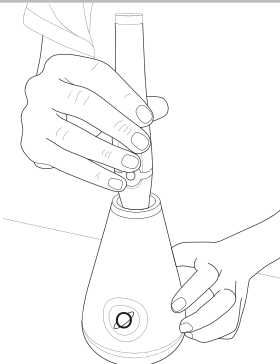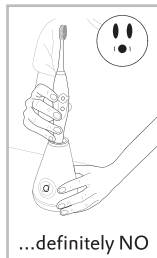


YES!



YES!

# Attention Please!

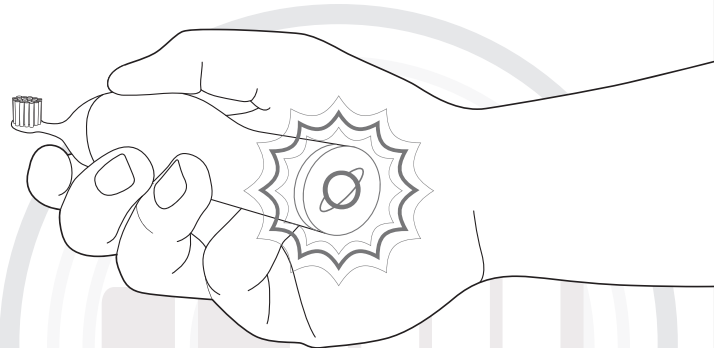


## Low Battery

The low battery alert is a beacon of distress for your Aura Clean Brush. The LED alert only triggers and pulses when the battery is low. Think of it as "Please feed me". If this occurs place the brush straight back into the Base Station. Give it a few moments and everything will be ok.

If the brush continues to pulsate while in the Base Station, please check to see that the power is connected.

# Aura Clean Dashboard

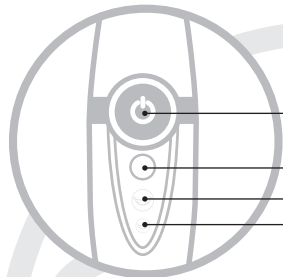

The Aura Clean Brush Dashboard puts you in control!


Power Button

Mode Selection Button

Super Clean Mode Indicator

Sensitive Clean Mode Indicator

## What's happening on the Aura Clean Dashboard?

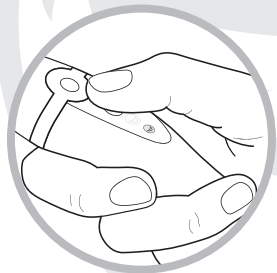

Press once for **Sensitive Clean Mode** and the small smiley will light up. Nice!




Press again to switch back to **Super Clean Mode** and the large smiley will light up. Easy!

# Activating the Brush


1

Remove the brush from the base

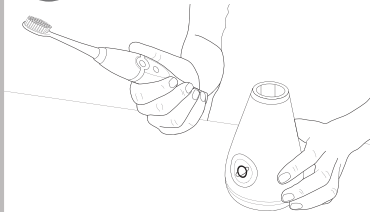

2


Select your mode, add your favourite toothpaste, water and press the power button

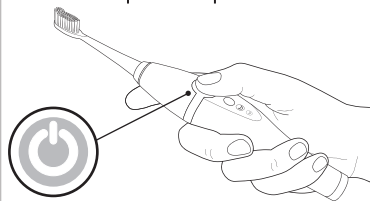

# Brushing Modes

Your Aura Clean Brush Handle has two cleaning modes and automatically starts in our recommended **Super Clean mode**.

**Sensitive Clean Mode** can be accessed at a touch of a button. **Sensitive Clean Mode** offers a gentle, yet thorough cleaning experience for sensitive gums and teeth. The clever brush even remembers the last mode set, so no need to change it.

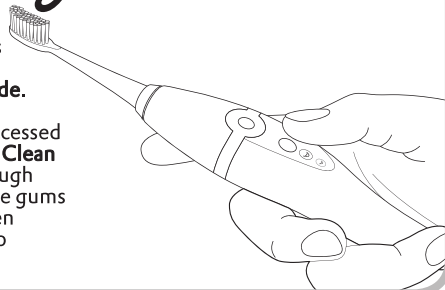

# Brushing Technique

Wet the bristles of your Aura Clean Brush and apply a small amount of your favorite toothpaste.


Now gently place the toothbrush bristles against the teeth at a slight angle towards the gum line. Press the power button and fire up the engine of the Brush. Sounds great doesn't it?

**Trust your Brush:** It will do all the brushing for you. Sonic technology is highly regarded and established as the next step evolution from the regular manual toothbrush. Its vibrational nature delivers thousands upon thousands of tiny brush strokes to safely dislodge debris, hygienically clean and efficiently polish. Think of it as your own dedicated army of 1400 personal assistants scrubbing away all the dirt from a ceramic wall. They get big job satisfaction from a shiny wall. Your teeth are that wall. The brush is happy to clean, so do not press hard or scrub.

**Hint:** Hold the brush head in place for a few seconds before sweeping on to the next tooth. It's saying hello and then goodbye. Allow the brush to work the gums as well as the teeth. Try to imagine yourself as an orchestra conductor, guiding cleanliness in your mouth.

The Aura Clean Brush delivers the dentist recommended brushing time of two minutes. After the two-minutes of brushing, the device stops automatically. At every 30 second interval the Brush will alert you to change brushing areas with a brief pause.

We've created a fun map of your mouth to demonstrate how to use the intervals. Please study it. You could say the teeth are tiny symmetrical islands in an ocean but we won't as they are just teeth.

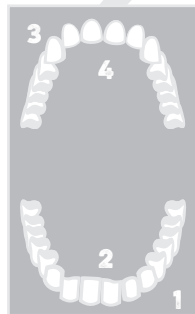

## The trusted TAO brushing technique

Follow our master guide for exceptional brushing. We've split the mouth into quadrants. Start brushing your outside bottom teeth (1) using the interval to move to your inside bottom teeth (2). After the pause, continue brushing the outside top teeth (3) until the last pause signifies the home stretch along the inside top teeth. You can actually mix it up, like toothbrush jazz, but don't miss a quadrant or play a quadrant twice. That's not the right cleaning rhythm.

After you've finished, feel free to conduct more additional cleaning if you need to remove stains or like the feeling. Just go easy and don't press.

We also recommend that you don't use any abrasive toothpaste that may hurt the tooth enamel.

Your Aura Clean Brush is safe to use on:

- Braces (brush heads wear out sooner when used on braces)
- Dental restorations (fillings, crowns, veneers)


# Brush Maintenance

## Keep Clean & Carry On

If you are using a filthy brush in your mouth, it isn't living up to its intended use.

Please follow these easy do maintenance routines to ensure your Aura Clean System is the very best condition to serve you rightly.

Keeping to these routines will also ensure your system is super happy.

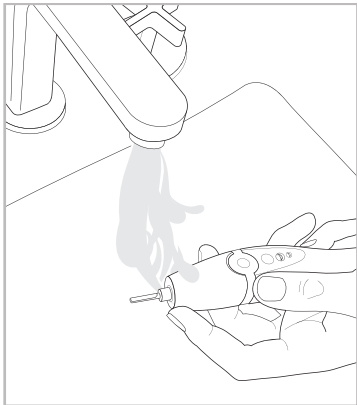

## Essential Cleanliness Maintenance Routines Keep the Aura Clean Brush Clean

Always rinse the Brush Head and bristles after each use.



We recommend that at least once a week you remove the brush head and rinse the metal shaft area with warm water.

Be careful not to damage the seal on the metal shaft with sharp objects when you clean.

# Base Station Cleaning

## Essential Cleanliness Maintenance Routines For The Base Station

We took this in consideration when we designed the smooth curved interior shape of the Base Station Core. We do recommend cleaning once a week to ensure tip top condition. Use a soft cloth. Wrap the cloth around your finger like so (see diagram no. 1). Unplug the Base Station. Carefully run your finger around the bottom (see no. 2) and do the same around the top-docking hole (see no. 3).

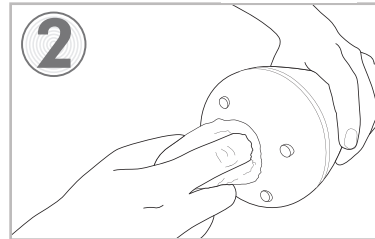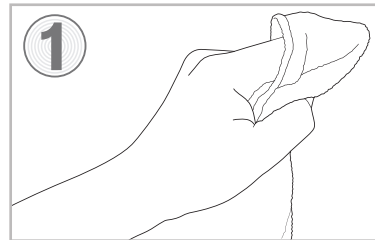

# Frequently Asked ...

Our team of Aura Clean experts have compiled a list of the questions most frequently asked about the Aura Clean System.

If you cannot find the answer to your question, please contact the TAO Clean Care Center in your country.

## **Why doesn't the Aura Clean Brush toothbrush work?**

You may need to recharge the Brush. Please remember the Brush was issued with enough energy for one brush only. Please charge the Brush for at least 24 hours before additional use.

Make sure the adaptor is plugged into a live electrical outlet and not into an outlet controlled by a light switch.

Is the brush head fixed on properly? You may need to reattach or replace the brush head (see page 21).

## **Why doesn't the UV sanitizer light go on?**

Did you place the Brush into the Base Station? Also make sure that the Brush is properly located within the docking hole.

Please check the wiring. Maybe you need to plug the Base Station into a live outlet.

Also, the Base Station needs to be on a smooth counter for the light to turn on.

## **Why does the toothbrush seem less powerful?**

You may need to charge the Brush.

You may need to replace the brush head. Make sure you change your brush head every 12 weeks for best results.

Perhaps it inadvertently went into sensitive mode.

## **Why is there a brief pause in the bristle motion while I'm brushing?**

Don't worry. This is the Brush reminding you to move to a new section of your mouth. These little reminders ensure you to spend equal time brushing all areas of the mouth (see page 18 for more info on technique).

## **The brush head rattles against my teeth. What can I do about this?**

Are you angling the brush correctly? To avoid this, keep the brush head bristles placed on the gum line at a slight angle. Perhaps you're pressing the hard plastic portion of the brush against your teeth. Don't do that, it's unpleasant at best.

## **Why do I feel a tickling sensation when I use the Aura Clean Brush?**

That's the feeling of being cleaned. You may experience a slight tickling/tingling sensation when you use Brush for the first time. Don't worry you'll get used to it. Just think how clean and shiny your teeth will be after a week.

## **Do I need a special charger if I travel with my Aura Clean System?**

The Aura Clean Base can be used at common worldwide voltages between 100-240V AC, at either 50 or 60Hz. All you need is a simple plug adapter. Also the Brush battery is designed to last up to 2 weeks between charges.

## **Does my Aura Clean System tell the future with its pulsing LEDs?**

Well, er no. We haven't quite got there with that particular technology. However, if you make sure your System is well maintained and used properly, it will keep you clean well into the future.

# Warranty

## Two year warranty

The Aura Clean System is warranted against defects in materials of workmanship for two years for the original purchaser from the date of purchase. It does not cover damages or wear resulting from the accident, misuse, abuse, commercial use, or unauthorized adjustment and/or repair. Please note: It does not cover wear resulting from using the brush heads.

There are no expressed warranties except as listed above. This warranty gives you specific legal rights, and you may have other rights which vary from state to state.

Should your Aura Clean System require service (or replacement at our option) while under warranty, do not return to the retailer. Please contact our customer service. Our friendly staff are trained to help you. If your Aura Clean System needs to be returned we will arrange a return label.

For assistance, please call +1-218-382-5326  
or email: [support@taoclean.com](mailto:support@taoclean.com)



Retailer stamp





## CONTACT

---

Further information

For assistance

call or text: +1-218-382-5326

email: [info@taoclean.com](mailto:info@taoclean.com)

---

T A © clean © 2015

All rights reserved. No part of this book may be reproduced, stored in a retrieval system or transmitted in any form or by any means, electronic, mechanical or otherwise, without written permission of the copyright owners.

