IC0907

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MASSAGE ACCESSORY  
ACCESSOIRE DE MASSAGE

Thank you for purchasing this product. Please read this instruction manual carefully before operating it.

Please pay attention to the important safety warnings and keep this manual for future use. Note: EBI reserves the right to revise the design and description of this product without notice.

## CONTENTS

Safety and maintenance .....	2
Features and benefits .....	3
Components .....	4
Operation.....	5
Applications.....	6
Reflexology Chart with corresponding body parts.....	7
Product specifications.....	8
Warranty .....	8
<b>FRANÇAIS .....</b>	<b>9</b>

2

## 1. SAFETY WARNINGS

- This product is not recommended for children.
- Disconnect the massage unit before cleaning
- Please operate this product according to the instructions in this manual.
- This product is designed for indoor use only.
- Recommended use time is 30 minutes or less.
- Please avoid use if the massage unit is damaged and call your iComfort technician or dealer.
- Do not use wet hands to pull the plug.
- Do not spill liquids onto this product to avoid electric shock or cause damage to this product.
- In the event of a power failure disconnect the massager

## 2. ENVIRONMENT

The iComfort massager should be used only in well ventilated, normally dry and dust free environments.

## 3. RESTRICTIONS

Stop immediately the use of this massage unit if you feel unwell and consult your doctor.

## 4. MAINTENANCE

- This product should only be maintained by iComfort certified technicians otherwise, the warranty will be void.
- Make sure the wall outlet is in good working order before use.
- If this massage unit will not be used for an extended period of time, store in a dry and dustless environment.
- Avoid product exposure to high temperatures and direct sunlight.
- Clean this product with a dry cloth. Do not use thinner, benzene or alcohol.
- We recommend to use the massager intermittently and not for long time periods.

## 5. TROUBLESHOOTING

- It is normal for the motor to make noise while in use.
- If the controller doesn't work normally, please check that the plug and socket are properly connected and the switch is turned on.
- When the programmed working time is over, the massager will power-off automatically.
- After several use, the massager might power-off because of overheat. It can be used again after half an hour of cooling time.

## FEATURES AND BENEFITS

The new iComfort ic 0907 foot massager foot and body provides strong, fast, deep-down relief. Just rest your feet, calves, thighs or lower back against the iComfort ic0907 foot and body massager and relax. Caress and soothe tired aching feet with a powerful vibrating action or use its multiple reflexology pressure points and infrared heating to stimulate blood circulation even more. The foot massage of your life is just a button-push away with its included wireless remote control. This unit is contoured to find known acupressure nodes and bring relief after the strains of a long day.

### Features:

- Automatic 15 minute timer
- Wireless remote control
- 8 Infrared heating massage levels
- 15 speed adjustment levels
- LED display


4

## COMPONENTS

### CONTROL PANEL
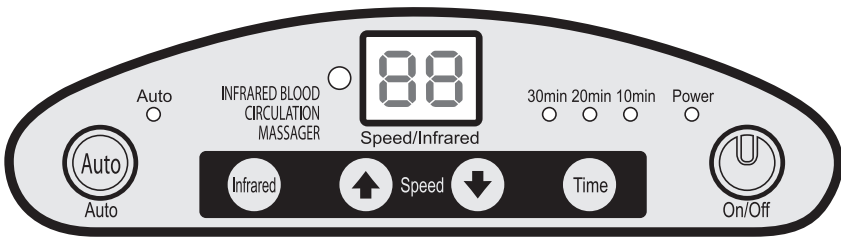
### REMOTE CONTROL

## OPERATION


- Securely position the iComfort massager in a stable place and connect it to an AC wall outlet.
- Press the “on” switch. The power light will go on and the massage session will start. On the panel, the timing light will come on, setting the timer at 10 minutes. Simultaneously, the LCD display will show 10, indicating 10 minutes. Pressing a second time the ‘time’ key will set the timer at 20 minutes. Pressing a third time the ‘time’ key will set the timer at 30 minutes. When the timer is about to stop, the indicator will flash. Press the “▲” key to increase the massage speed, and the “▼” key to decrease it. You can choose between 0-15 to select the speed level of your choice.
- There are eight levels of infrared heat to choose from: H0-H8. H0 will turn the infrared heating off.
- To select an automatic pre-programmed massage mode, press the ‘automatic’ key and the iComfort massager will massage according to different sets of massage speeds. Pressing the “speed” key will automatically display the previous preset data.
- The programmed massage length, speed and heat setting should vary according to your health. We recommend to start using the iComfort massager slowly and progressively as you get use to the unit more and more.



## APPLICATIONS

In addition to the massage of the feet, the iComfort infrared foot massager may also be used to massage the upper and lower back, your thighs, calves, and hands.

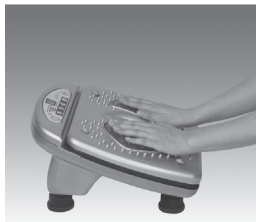**BACK****WAIST****ABDOMEN****FEET****PALM****CALVES**

### REFLEXOLOGY CHART WITH CORRESPONDING BODY PARTS

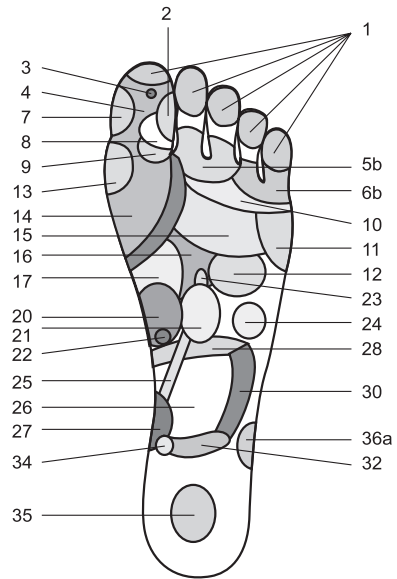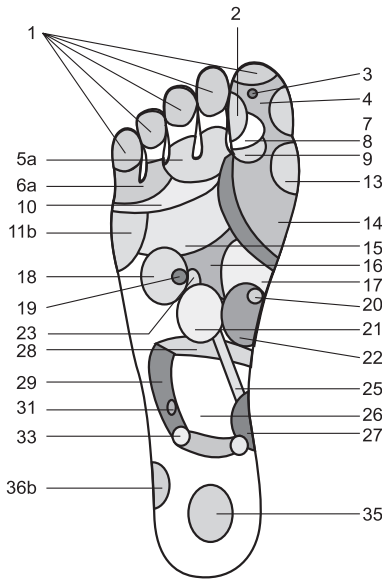- 1. Forehead
- 2. Sun point
- 3. Pituitary gland
- 4. Cerebrum
- 5a. Left eye
- 5b. Right eye
- 6a. Left ear
- 6b. Right ear
- 7. Nose
- 8. Cerebellum
- 9. Neck
- 10. Muscle at the back of neck
- 11. Shoulders
- 12. Heart
- 13. Accessory thyroid glands
- 14. Thyroid glands
- 15. Bronchial tubes
- 16. Solar plexus
- 17. Stomach
- 18. Liver

- 19. Gall bladder
- 20. Pancreas
- 21. Kidneys
- 22. Duodenum
- 23. Adrenal glands
- 24. Spleen
- 25. Ureters
- 26. Small intestine
- 27. Bladder
- 28. Transverse colon
- 29. Ascending colon
- 30. Descending colon
- 31. Ileum entrance
- 32. Rectum
- 33. Appendix
- 34. Anus
- 35. Ovaries
- 36a. Left knee
- 36b. Right knee

## PRODUCT SPECIFICATIONS

Rated Voltage: AC110V 60Hz  
 Rated power: 65 W

## WARRANTY

### LIMITED ONE (1) YEAR WARRANTY

This product is guaranteed against all manufactured faults and defects for a period of one (1) year from the date of purchase.

For repairs or product defects, please call our customer service department at : 1-866-967-7333. Please have your proof of purchase on hand. If we are not able to replace your defective part, we will send you a service technician for repairs.

Important notes: Any damages caused by an accident or an abusive usage is NOT covered under warranty. We also decline all responsibility concerning indirect or consecutive damages. In certain Canadian provinces, it is not permitted to limit indirect or consecutive damages. Therefore it is possible these limitations are not applicable to you. This warranty gives you specific legal rights. You may also have other legal rights which vary from one province to another.

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Pictures are for reference only, colors and models may vary without notice. Made in China. iComfort™ is a trademark of EBI Inc. Photographies pour référence seulement. Les couleurs et modèles sont sujet à varier sans préavis. Fabriqué en Chine. iComfort™ est une marque de commerce de EBI inc.